

A close-up, profile view of a person's face, focusing on the mouth and teeth. The person is smiling, showing a row of bright, white teeth. The lips are a natural pink color. The background is a solid light blue.

Dental Products & Services Catalog

Our Mission

Our mission is to provide our customer, the dental professional, with products and services that make them more successful in the eyes of their customer, the dental patient.

At DenMat, we develop, manufacture, and directly distribute dental equipment, materials, instruments, lab services, and provide education to help dental professionals deliver beautiful smiles to patients worldwide... because we believe everyone deserves a healthy, beautiful smile.

Our Commitment to You

Our customers come first. Customer “service” is not enough for us. Our paramount goal is to ensure that doing business with DenMat is not only easy, but a pleasure.




DenMat's state-of-the-art global headquarters is located in Lompoc, CA, in the heart of California's Central Coast. Our 100,000 square-foot facility houses most of our North American operations, including manufacturing, lab, sales, marketing, warehousing, distribution, administrative services and executive offices.




State-of-the-Art Learning Center



State-of-the-Art Production

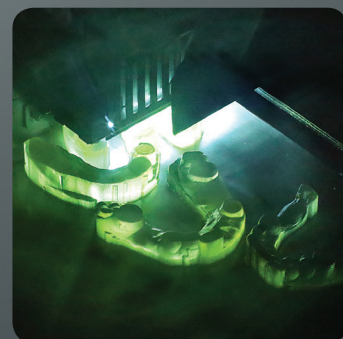


State-of-the-Art Lab Operations

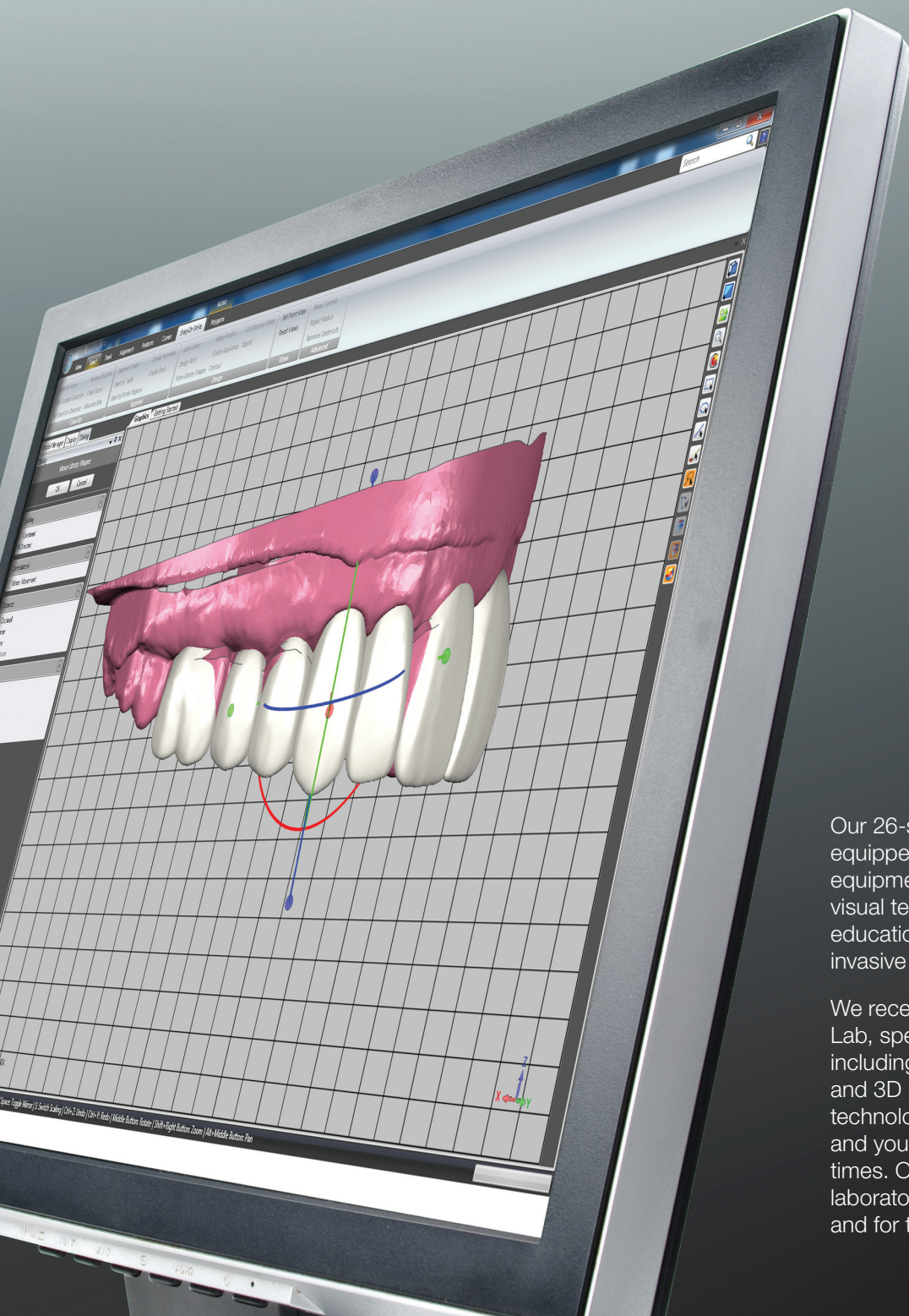

Our 26-seat Dental Learning Center is equipped with the latest in dental operator equipment, along with high-definition audio/visual technology for hands-on and live-patient education on Lumineers and other minimally invasive esthetic dental programs.

We recently invested millions into the DenMat Lab, specifically in technology and equipment including scan and design software, 5-axis mills and 3D printers. This investment in CAD/CAM technology means consistent quality for you and your patients, along with faster turnaround times. Our goal is to be your first-choice dental laboratory for minimally-invasive esthetics today and for the future.

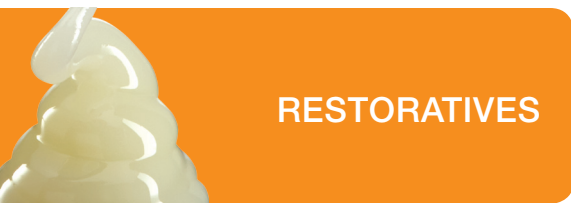

RESTORATIVES

Achieve outstanding esthetic results using DenMat's award-winning restorative products. Whether a cement, adhesive, or composite, dentists have been relying on DenMat for over 40 years as the source for superior restorative materials.




IMPRESSIONS

Choose DenMat materials for reliable, consistent impressions, every time. From simple cases to complex full-mouth restorations, our award-winning VPS materials and impression trays deliver the quality and performance you depend on.



WHITENING

DenMat's comprehensive line of both in-office and take-home whitening products will help you build your practice while delivering extraordinary results you and your patients will love.

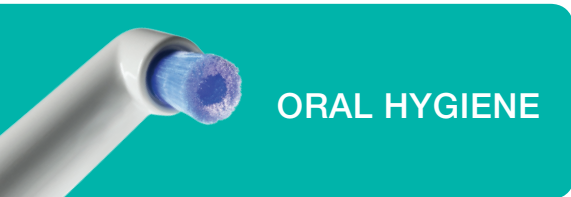

ORAL HYGIENE

DenMat's oral hygiene line features the Rotadent ProCare™ power toothbrush, Pro-Select Platinum® Piezo ultrasonic scaler, ATRIDOX® locally-applied periodontal antibiotic and Pro-DenRx® fluoride gels, rinses and varnish. A two-level Soft Tissue Management® (STM®) CE offering rounds out this portfolio that is dedicated to the detection, prevention and treatment of periodontal disease.



INSTRUMENTS

DenMat acquired G. Hartzell & Son in 2014. These industry-leading precision products have long been preferred by clinicians for perio, hygiene, restorative, ortho, endo, and surgical dental procedures. Simply put, their quality is unparalleled.



MAGNIFICATION & ILLUMINATION

PeriOptix® has been manufacturing top-of-the-line loupe and headlight systems since the late 90s, a tradition that continues today after being acquired by DenMat in 2013. All our innovative loupes and headlights are manufactured in California and are made of cutting-edge components with strict quality standards.



SMALL EQUIPMENT

DenMat offers a diverse line of small equipment, including soft-tissue lasers, curing lights, and oral lesion detection systems. Carefully engineered and rigorously tested, our products have become must-haves in dental offices all over the world.





DENMAT LAB

Focused on esthetics, the DenMat Lab produces a full range of restorations for just about any need. From simple, monochromatic smile improvement to lifelike AACD-quality work, we deliver esthetics and function with every case.

Table of Contents

Restoratives	4	Instruments	40	Magnification & Illumination	50
Core Build-Up	6	Hygiene Instruments	42	Firefly Cordless Headlight System	52
Core Paste.....	6	Scaling & Root Planing Kit.....	42	LED Headlamps	53
Core-Post.....	7	Rigid Gracey Curettes.....	43	Firefly Cordless Headlight System.....	53
Temporary	8	Barnhart Universal Curette.....	43	MicroLine Mini LED Headlight.....	53
Perfectemp10.....	8	Columbia Universal Curette.....	43	Solaris Portable LED Headlight.....	53
Perfectemp II.....	9	UC/Younger Good Universal Curette.....	43	PeriOptix Through-the-Lens Loupes	54
Adhesives	10	Loma Linda Curette.....	43	2.3x MicroLine TTL Loupes.....	54
Tenure4G.....	10	Implant Curettes.....	43	2.7x TTL Loupes.....	54
Tenure MPB.....	11	Sickle Scalars.....	43	3.1x TTL Loupes.....	54
Tenure Uni-Bond.....	11	Explorers.....	43	PeriOptix Flip-Up Loupe Optics	55
Tenure Quik.....	11	Periodontal Probes.....	43	Panoramic Flip-Up Loupes.....	55
Gloss-N-Seal.....	11	Cone Socket System.....		Legacy Sport Flip-Up Loupes.....	55
Cements	12	Cone Socket Instruments.....	43	MicroLine Flip-Up Loupes.....	55
Infinity.....	12	Universal Curettes.....	43	Frames	56
Geristore.....	14	Sickle Scalars.....	43	PeriVista Frame by Adidas.....	56
UltraBond Clear.....	15	Rigid Gracey Curettes.....	43	Legacy Sport Frames.....	56
Ultra-Bond Plus.....	15	Restorative Instruments	44	PeriOptix Flip-Up Loupes/Frame Systems	57
Composites	16	Composite Kit.....	44	Air-X Flip-Up Loupes.....	57
Nuance Universal.....	16	Spoon Excavators.....	45	Hogies Advanced Optical Platform.....	57
Nuance FLOW.....	17	Excavators.....	45	Small Equipment	58
Virtuoso Universal.....	17	Explorers.....	45	NV Microlaser.....	60
Virtuoso Flowable.....	17	Cord Packers.....	45	SOL.....	62
Impressions	18	Carver Cord Packer.....	45	Sapphire.....	63
SplashMax.....	20	Carver.....	45	Flashlite Magna 4.0.....	63
Splash!.....	22	Hollenback Carver.....	45	Flashlite 2.0.....	63
Precision.....	23	Interproximal Carver.....	45	Lab	64
Bite Registration.....	24	Discoïd-Cleoid.....	45	Lumineers.....	66
The Gripper.....	25	Ball Burnishers.....	45	LumiZir.....	66
Perfectray.....	25	Composite Instruments.....	45	Snap-On Smile.....	66
Whitening	26	Crown Removers.....	45	A.R.T. (Additive-Reductive Template).....	67
LumiSmile White Take-Home.....	28	Articulating Paper Forceps.....	45	Signature Service.....	67
LumiSmile White In-Office.....	29	Surgical Instruments	46	Lumineers & Snap-On Smile.....	
Lumibrite Take-Home.....	29	Implantology Starter Kit.....	46	Kits/Replacements.....	67
Lumibrite Chairside.....	29	Periotomes.....	47	Professional Dental Photography Kit.....	67
Oral Hygiene	30	Periosteals.....	47	Lab Product – Finishing & Polishing	68
Power Toothbrushes	32	Surgical Curettes.....	47	LumiGrip.....	68
Rotadent ProCare.....	32	Needle Holders.....	47	Ceri-Saw & Ceri-Sander.....	69
Ultrasonic Scalars	34	Curved Scissors.....	47	Lumineers Finishing Kit.....	69
Pro-Select Platinum (PSP).....	34	Perio Knives.....	47	Matrix Strips.....	69
Preventatives	36	Elevators.....	47	Practice Enhancement & Education	70
Pro-DenRx 2.0%.....		Root Tip Picks.....	47	LSDP (Lumineers Smile Discovery Program)	72
Neutral Sodium Fluoride Rinse.....	36	Luxating Elevators.....	47	LumiSmile.....	72
Pro-DenRx 0.2%.....		Scalpel Handles.....	47	National Advertising.....	73
Neutral Sodium Fluoride Rinse.....	36	Extraction Forceps.....	47	LSDP Member's Website.....	73
Pro-DenRx 0.63%.....		Deep-Gripping Extraction Forceps.....	47	LumiDays.....	73
Stannous Fluoride Rinse.....	36	Endodontic Instruments	48	Continuing Education	74
Pro-DenRx 0.044% Daily APF Rinse.....	36	Endodontic Starter Kit.....	48	Soft Tissue Management (STM).....	74
Pro-DenRx Fluoride Varnish.....	37	Endo Explorers.....	49	Mastering Minimally-Invasive Smile Design.....	75
Pro-DenRx Fluoride Gels & Pastes.....	37	Endo Excavators.....	49	E-Commerce	76
Pro-DenRx Chlorhexidine Gluconate.....	37	Endo Condensers.....	49	Terms & Conditions	78
Oral Hygiene Essentials	38	Periosteal.....	49	Contact Us	80
ATRIDOX.....	38	Dovgan MTA Carrier.....	49		
ViziLite TBlue.....	39	Terauchi Gutta-Percha Removal Instruments.....	49		
RotaPoint.....	39	Dovgan MTA Pluggers.....	49		
Pro-DenRx Probes.....	39	Perf Repair Instrument.....	49		
Pro-Ties.....	39	Huey Forceps (Diamond Coated).....	49		
		Peets Forceps.....	49		
		Micro Mirrors.....	49		
		Simple Stem Handle.....	49		

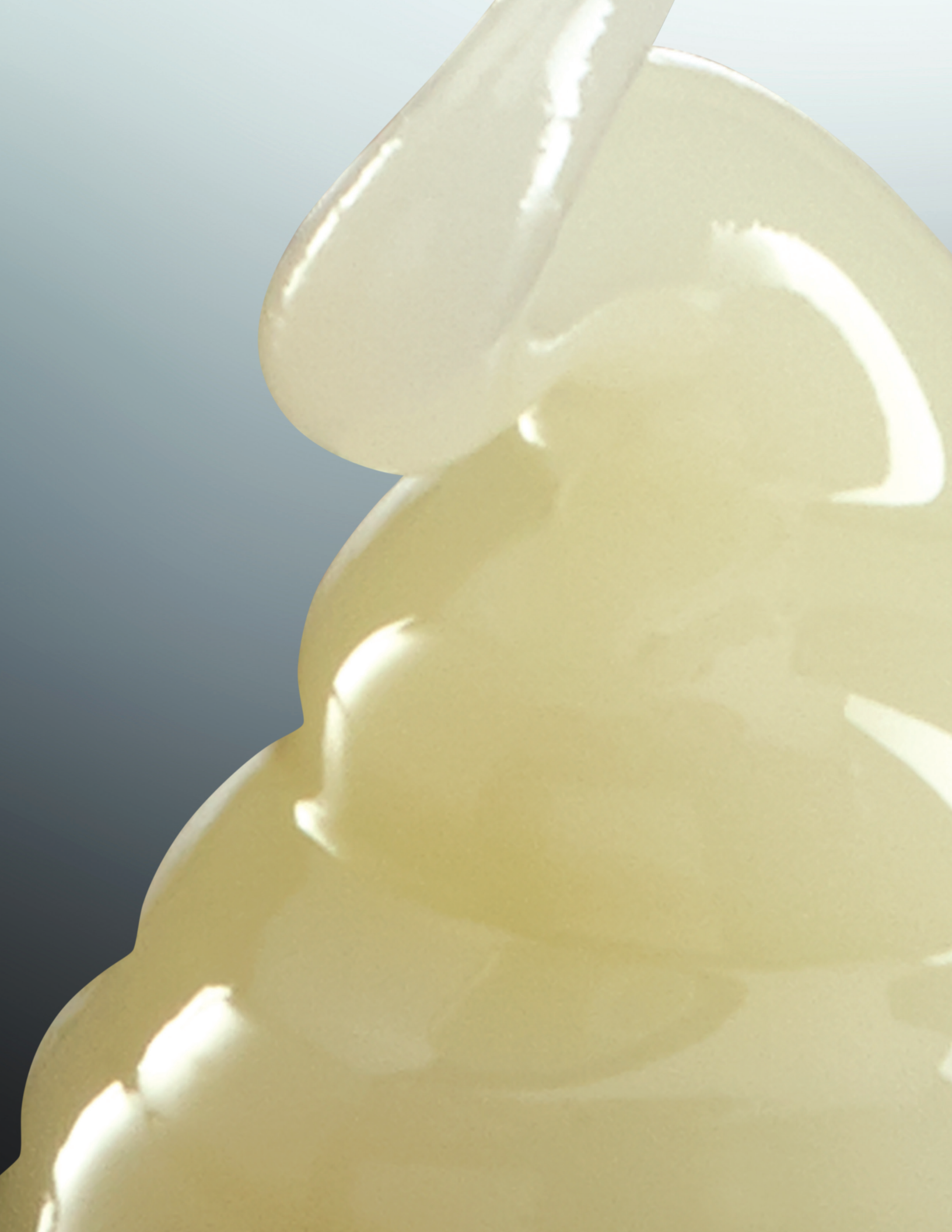

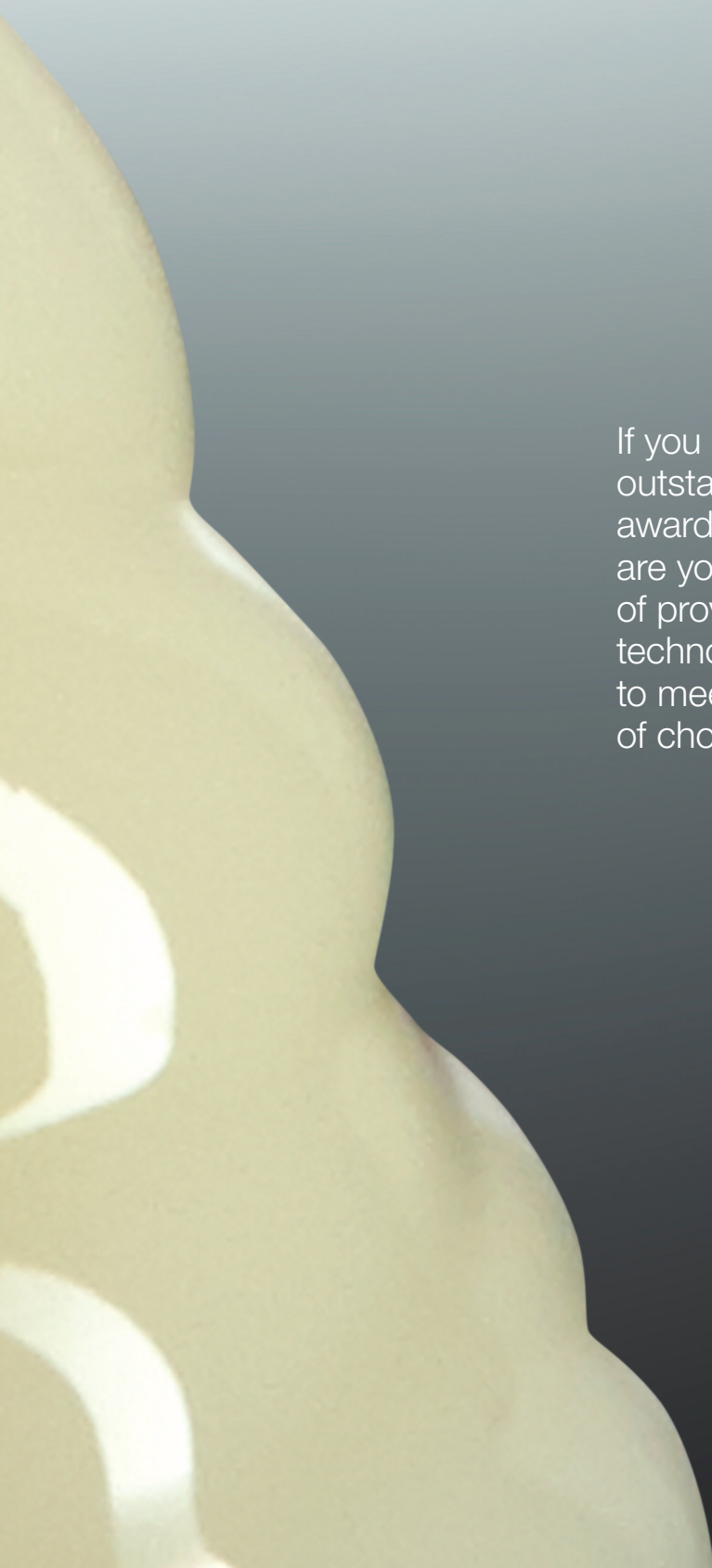An abstract graphic on the left side of the page shows a cross-section of a tooth. The crown is a light yellow color, and the root is a darker yellow. The background is a gradient from light blue at the top to dark grey at the bottom.

If you are seeking a way to achieve outstanding esthetic results, DenMat®'s award-winning restorative products are your solution. Our extensive line of proven materials and advanced technology offer clinicians solutions to meet their style and technique of choice.

CORE BUILD-UP



Core Paste

The Proven Time Saver for Core Build-Ups and Posterior Restorations

For over 30 years, Core Paste has helped thousands of dentists perform reliable core build-ups with complete confidence.

- Dual-cure formula
- Stackable viscosity and excellent handling characteristics
- Cuts like dentin
- Radiopaque




Before



After

Core Paste's stackable viscosity allows build-up without a matrix or core form!

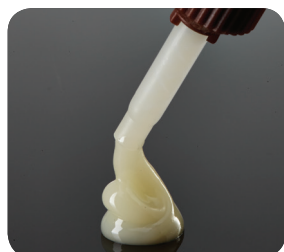

Product	Item #
Core Paste XP	
Core Paste XP Enamel with Fluoride Intro Kit*	030657110
Core Paste XP Dual-Cure Enamel with Fluoride*	030652110
Core Paste XP Dual-Cure Enamel	030653110
Core Paste XP Dual-Cure White	030654110
Core Paste XP Dual-Cure White with Fluoride*	030655110
<i>Contents: 4 – 10g syringes (Intro Kit contains 2 syringes only), auto-mixing tips (40pc), clear intraoral tips (40pc), intraoral tips 18G (6pc), Tenure A & B bonding agent, Etch+N'Seal, light cure test rings (3pc)</i>	

"I've been using Core Paste in my practice for many years. The stackable viscosity provides excellent handling properties. Core Paste is easy to use and consistently results in strong and effective restorations."

— Ian Shuman, DDS, Pasadena, MD

Product	Item #
Core Paste Syringeable	
Core Paste White with Fluoride Self-Cure*	030625100
Core Paste White Self-Cure	030635100
Core Paste White with Fluoride Self-Cure & Gun*	030625500
Core Paste White Self-Cure & Gun	030635500
Core Paste Enamel Dual-Cure	030645100
Core Paste Enamel Dual-Cure & Gun	030645500
<i>Contents: 1 – 50g cartridge, auto-mixing tips (20pc), intraoral tips (20pc), light cure test rings (3pc), mixing pad, mixing tips, instructions</i>	

Product	Item #
Core Paste Jars	
Core Paste White Q/C with Fluoride Kit* – 50g	030612525
Core Paste White with Fluoride Kit* – 50g	030622100
Core Paste Enamel Kit – 50g	030616000
Core Paste Enamel with Fluoride Kit* – 50g	030616025
Core Paste Enamel Q/C Kit – 50g	030616325
Core Paste Blue S/C Kit – 50g	030617325
Core Paste White Kit – 50g	6100
<i>Contents: 2 – 25g jars Core Paste (A & B), centrix tip/plunger orange/green (25pc), mixing pad, mixing sticks, all-surface bonding chart, instructions</i>	



Core-Post

Ideal for Core Build-Ups in Single-Tooth Restorations

Flexing like natural tooth structure, Core-Post dramatically minimizes tooth stress and fractures that occur with standard metal posts. Core-Post is available in both carbon fiber for increased posterior strength and tooth-colored glass fiber for anterior esthetics.

- Easy to cut and use
- Radiopaque
- Available in carbon fiber and glass fiber
- Five diameters: 1.0, 1.2, 1.4, 1.7 and 2.0mm

Product	Item #
Ancillary Core Build-Up Items	
Centrix Tip/Plunger Orange/Green (25pc)	030057001
Auto-Mixing Tips LG ORF Blue (20pc)	030594121
Intraoral Tips 25ml Yellow (20pc)	030594101

Product	Item #
Ancillary Core Build-Up Items	
Syringe Mixing Tips Brown/Clear (10pc)	030777615
Clear Intraoral Tip 10G (10pc)	030653105

Product	Item #
Core-Post Value Kit	
Core-Post Value Kit	033637000
<i>Contents: 5 – 1.0mm, 1.2mm, 1.4mm, 1.7mm, 2.0mm glass fiber posts, 5 – 1.0mm, 1.2mm, 1.4mm, 1.7mm, 2.0mm carbon fiber posts, all-surface bonding chart, instructions</i>	
Core-Post Intro Kit	
Core-Post Intro Kit	033625000
<i>Contents: 1 – 1.0mm, 1.2mm glass fiber posts, 1 – 1.4mm, 1.7mm carbon fiber posts, instructions</i>	
Core-Post 10-Piece Kits	
Core-Post Glass 1.0mm Kit	033638010
Core-Post Glass 1.2mm Kit	033639010
Core-Post Glass 1.4mm Kit	033640010
Core-Post Glass 1.7mm Kit	033641010
Core-Post Glass 2.0mm Kit	033642010
Core-Post Carbon 1.0mm Kit	033643010
Core-Post Carbon 1.2mm Kit	033644010
Core-Post Carbon 1.4mm Kit	033645010
Core-Post Carbon 1.7mm Kit	033646010
Core-Post Carbon 2.0mm Kit	033647010

dm
DenMat®

TEMPORARY



Perfectemp10

Beautiful Temps. Beautiful Price.

Perfectemp10 is a premium multi-functional acrylic composite temporary crown and bridge material that delivers the strength that patients need at a competitive price that dentists desire.


- 10:1 automix formulation
- Superior strength
- Natural-looking luster
- Cartridge or syringe
- Five Vita shades

**perfect
temp¹⁰**


Temporary Crown & Bridge Material



Take an impression of treatment area

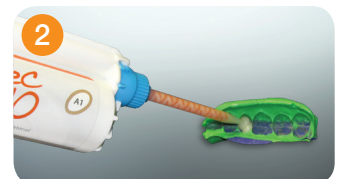

Extrude Perfectemp10 into impression

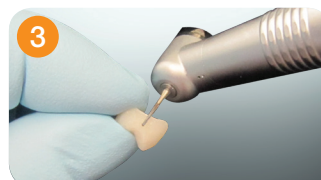

Trim, contour and finish temporary



Cement temporary in place



Product	Item #
Perfectemp10 Syringe – 10ml	
A1	CR1218
A2	CR1220
A3.5	CR1219
<i>Contents: 1 – 10ml syringe, mixing tips, instructions</i>	
Perfectemp10 Cartridge – 50ml	
A1	SPD1305
A2	SPD1306
A3.5	SPD1307
B1	SPD1308
BL-2	SPD1309
<i>Contents: 1 – 50ml cartridge, mixing tips, instructions</i>	

“Perfectemp10 requires no polishing or glazing to achieve a natural surface texture and dentition integration. It saves me time and creates exceptionally strong and beautiful temporaries.”

— Simona Cuevas, DDS, San Antonio, TX

Perfectemp II Strong and Esthetic Temporaries at a Practice-Friendly Price

Perfectemp II is a reliable, self-cure, bis-acrylic provisional crown and bridge material that produces durable and esthetic temporaries.

- 6-minute set time; 2.5-minute removal time
- Easy trimming and polishing for a natural luster
- Five Vita shades for a perfect match
- Automix
- Bis-acrylic




Product	Item #
Perfectemp II Syringe – 5ml	
A1	CR1039
A2	CR1043
A3.5	CR1042
B1	CR1041
BL-2	CR1040
<i>Contents: 2 – 5ml syringes, mixing tips, instructions</i>	
Perfectemp II Cartridge – 50ml	
A1	SPD1300
A2	SPD1301
A3.5	SPD1302
B1	SPD1303
BL-2	SPD1320
<i>Contents: 1 – 50ml cartridge, mixing tips, instructions</i>	



Dispensing
Gun

Product	Item #
Mixing Tips	
Cartridge Mixing Tips (25pc)	SPD1304
Syringe Mixing Tips (25pc)	CR1046

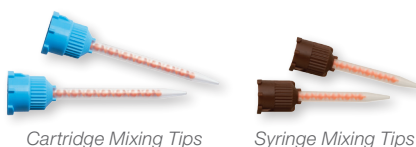

Cartridge Mixing Tips

Syringe Mixing Tips

Product	Item #
Dispensing Guns	
Dispensing Gun for Cartridges 4:1 (Perfectemp II)	SPD1178
Dispensing Gun for Cartridges 10:1 (Perfectemp10)	034381150
Dispensing Gun Insert for Cartridges 10:1	040381150



Tenure4G

Dual-Solvent Bonding Agent

4Greater bond strength and sensitivity control, turn to Tenure4G, DenMat's next generation of the popular Tenure MPB bonding system. Tenure4G is a 4th generation multi-purpose, self-cure adhesive system for bonding any resin restorative to all intraoral surfaces. Tenure4G's formulation contains added ethanol to reduce patient sensitivity. Additionally, it has shown improved bond strength, to give you more confidence in your composite restorations.

- Ethanol and acetone for strength and sensitivity control
- Improved bond strength
- Self-cure
- Easy to use



tenure[®]4G



Tenure MPB Multi-Purpose Bonding System

Tenure Multi-Purpose Bonding System has high bond strength and excellent retention, along with significant reduction of microleakage and post-operative sensitivity, making it the bonding agent of choice for thousands of clinicians.

- Bonds to all surfaces
- Versatility
- Complete self-cure
- Easy flow
- Low film thickness



Tenure Uni-Bond® Self-Etching Bonding and Sealing System

Tenure Uni-Bond self-etching bonding system simplifies all direct bonding applications with anterior and posterior composites. It is ideal when cementing crowns and bridges with self- and dual-cure resin cements. Tenure Uni-Bond is compatible with self-cure, dual-cure and light-cure composite materials.

- No separate etching required
- High bond strength to all intraoral surfaces
- Less chair time
- Non-sensitivity formula
- Gloss-N-Seal seals and glosses final restoration



Tenure Quik® Single-Bottle System

Tenure Quik's formulation provides a surface-penetrating, low-viscosity bonding agent. Once applied to etched enamel and dentin, you can be assured that Tenure Quik with Fluoride will completely penetrate and seal the tooth interface.

- High bond strengths to all surfaces
- Single component saves time
- Simplified placement
- Compatible with self-cure, dual-cure and light-cure composite restorations



Gloss-N-Seal® Light-Cure Sealant

Gloss-N-Seal is a light-cure, low-viscosity, surface penetrating sealant with exceptional abrasion resistance. It penetrates and seals to reduce composite wear, marginal staining and breakdown.

- Sealing composite restorations
- Resurfacing old composite restorations
- Coating provisional for instant glaze finishes

Product	Item #
Tenure 4G	
Tenure 4G System	TEN4000
<i>Contents: 1 – 6ml A bottle, 1 – 6ml B bottle, 1 – 6ml Tenure S Bond Enhancer, 1 – 3ml Tenure S Filled Enhancer, 1 – 3ml Etch'N'Seal, syringe tips 18G (6pc), brush tips (60pc), dappen dish (2pc), brush-tip handle, all-surface bonding chart, instructions</i>	
Tenure MPB	
Tenure MPB System	031146000
<i>Contents: 1 – 6ml A bottle, 1 – 6ml B bottle, 1 – 6ml Tenure S Bond Enhancer, 1 – 3ml Tenure S Filled Enhancer, 1 – 3ml Etch'N'Seal, syringe tips 18G (6pc), brush tips (60pc), dappen dish (2pc), brush-tip handle, all-surface bonding chart, instructions</i>	
Tenure MPB A 6ml Refill Kit	030114510
Tenure MPB B 6ml Refill Kit	030114610
<i>Contents: 1 – 6ml bottle, instructions</i>	
Tenure Uni-Bond	
Tenure Uni-Bond Kit	030411900
<i>Contents: 1 – 6ml A bottle, 1 – 6ml B bottle, Gloss'N'Seal (3g jar), dappen dishes (2pc), micro brush (100pc), all-surface bonding chart, instructions</i>	
Tenure Uni-Bond A 6ml	030412901
Tenure Uni-Bond B 6ml	030412902
<i>Contents: 1 – 6ml bottle, instructions</i>	
Tenure Quik	
Tenure Quick with Fluoride Kit*	031145700
<i>Contents: 1 – 5ml bottle, 1 – 3ml etchant medium viscosity syringe, dappen dish (2pc), brush tips (25pc), brush-tip handle, pink syringe tips 18G (6pc), all-surface bonding chart, instructions</i>	
Tenure Quik with Fluoride Refill Kit*	031145300
<i>Contents: 3 – 6ml bottles, all-surface bonding chart, instructions</i>	
Tenure S	
Tenure S Value Kit	031145100
<i>Contents: 3 – 6ml bottles, instructions</i>	
Tenure S Bond Enhancer	032245101
<i>Contents: 1 – 6ml bottle, instructions</i>	
Tenure S Filled Bond Enhancer	031145110
<i>Contents: 1 – 3ml syringe, instructions</i>	
Dab-Eze Tenure S 3g Jar	031145201

Product	Item #
Gloss-N-Seal	
Gloss-N-Seal – 3g jar	030451010

Product	Item #
Ancillary Adhesive Items	
Syringe Tips Pink 18G/30 (100pc)	030007121
Syringe Tips Pink 18G/30 (20pc)	030007111



Infinity

All-Purpose Resin-Ionomer Cement

Infinity is the all-purpose resin-ionomer cement that is great for all zirconia and metal restorations. It is also ideal for porcelain-to-metal, all-ceramic, resin and ceramic inlays/onlays, cast metal and all types of posts.

- Cleans up in seconds
- Self-adhesive
- Dual-cure formula
- Low film thickness
- High bond strengths to all surfaces



Simply syringe Infinity directly into the crown and seat.



Removal of excess cement is fast and easy — no drilling or scraping required!

Technique for Bridge Cementation Using Infinity



Clean tooth surfaces thoroughly
Dispense Infinity into the restoration
Note: for inlays/onlays, syringe Infinity directly onto the tooth




Seat the restoration under occlusal forces



Briefly light cure the excess material extruded from the margins for 1 – 3 seconds until it reaches a gel consistency



Remove the gelled excess Infinity with an instrument

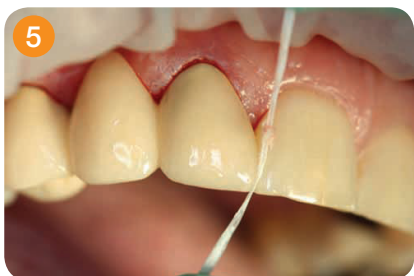

Use dental floss to remove excess material interproximally



Continue light curing process or allow the restoration to set for 2 – 3 minutes

Product	Item #
Infinity Intro Kit	
Infinity Intro Kit	039690110
Contents: 1 – 5g syringe, auto-mix dispensing tips (3pc), instructions	
Infinity Value Kits	
Infinity Syringeable Value Kit	039680100
Infinity SE Syringeable Value Kit	039681200
Contents: 3 – 10g syringes, auto-mix dispensing tips (60pc), instructions	
Infinity Syringeable Kits	
Infinity Syringeable Kit	039670100
Infinity SE Syringeable Kit	039680200
Contents: 1 – 10g syringe, auto-mix dispensing tips (20pc), instructions	
Infinity Jars	
Infinity Original A & B Jars	039620150
Contents: 1 – 8g A jar, 1 – 8g B jar, mixing pad, mixing sticks, instructions	

Product	Item #
Ancillary Cement Accessories	
Syringe Mixing Tips	030795525
8-Element (10pc)	

“Crown and bridge cement that is clinically proven to hold up beyond the test of time.”

— Dr. Louis Kaufmann, Chicago, IL





AMAZINGLY BIOCOMPATIBLE!

Tissue reattaches
to Geristore
as it would to
natural tooth
structure.

Geristore

Resin-Ionomer Restorative Material

Geristore is the self-adhesive, easy-to-use, resin-modified glass ionomer that you will find yourself using on your most difficult to your everyday restorations. Choose Geristore for a variety of clinical situations and many more:

- Base material/liner
- Pulp-capping agent
- Subgingival restorations
- Abfraction lesions/root resorption lesions
- Re-cementing crowns and bridges



Before: Abfraction lesions



After: Restorations successfully and esthetically completed with Geristore and veneered with a clear resin



Product	Item #
Geristore Syringes – Intro Kit	
Geristore Intro Kit	031476110
<i>Contents: 1 – 5g syringe shade A2, auto-mixing tips (6pc), intraoral tips (6pc), light cure test rings (3pc), instructions</i>	
Geristore Syringes – Value Kits	
Geristore Syringeable Value Custom Kit (choose 3 shades)	031477100
Geristore Syringeable Value Kit Shade A2	031458550
Geristore Syringeable Value Kit Shade A3	031457520
<i>Contents: 3 – 10g syringes, auto-mixing tips (60pc), intraoral tips (40pc), light cure test rings (3pc), instructions</i>	
Geristore Syringeable Kit	
Syringeable Kit Shade A2	031457500
Syringeable Kit Shade A3	031457720
Syringeable Kit Shade A3.5	031457730
<i>Contents: 1 – 10g syringe, auto-mix dispensing tips (20pc), intraoral tips (10pc), light cure test rings (3pc), instructions</i>	



Tenure & Ultra-Bond Plus
A winning combination for proven reliability and long-term performance!

UltraBond Clear Light-Cure Resin Cement

UltraBond Clear is a light-cure resin cement formulated for bonding Lumineers, as well as other porcelain or composite veneers. Featuring a viscosity developed for optimal handling when placing veneers, and highly transparent clear and natural shades, UltraBond Clear is sure to provide the final esthetics you and your patients will love.

- Easy-to-use, no-mix, light-cure application
- Excellent color stability

Ultra-Bond Plus Restorative Resin Cement

Ultra-Bond Plus is a multi-purpose, dual-cured, ADA Type II restorative, as well as a cement. Ultra-Bond was custom designed to match the physical properties of Cerinate Porcelain and tooth structure for long lasting restorations.

- Prevents marginal staining and pop-offs
- Creates long-lasting, durable restorations
- Ideal cementing/luting medium for crowns/veneers
- Protects against wash-out
- Exceptional ease of use



Clear Intraoral Tip

Pink Syringe Tip

Product	Item #
Ancillary Cement Accessories	
Syringe Mixing Tips Brown/Clear (10pc)	030777615
Clear Intraoral Tip 10G (10pc)	030653105
Syringe Tips Pink 18G/30 (100pc)	030007121

Product	Item #
Geristore Jars	
Geristore A Initiator – 10g	4506
Geristore A Initiator Slow Set – 10g	030450911
Geristore B Paste Shade A1 – 10g	034522101
Geristore B Paste Shade A2 – 10g	034523101
Geristore B Paste Shade A3 – 10g	034524101
Geristore B Paste Shade A3.5 – 10g	034526101

Product	Item #
UltraBond Clear Intro Kit	
UltraBond Clear Intro Kit	UB2000
<i>Contents: 1 – 3g syringe Clear Shade, 1 – 3g syringe Natural Shade, 1 – 3ml syringe Etch'N'Seal, 1 – 6ml bottle Tenure S Bond Enhancer, 1 – 6g bottle Porcelain Conditioner, 1 – 5ml bottle Cerinate Prime, syringe tips 18G (20pc), brush tips (60pc), dappen dish (2pc), brush-tip handle, instructions</i>	
UltraBond Clear Refill Kits	
UltraBond Clear Refill Kit – Clear Shade	UB2001
UltraBond Clear Refill Kit – Natural Shade	UB2002
<i>Contents: 3 – 3g syringes, syringe tips 18G (30pc), light curing test rings (3pc), instructions</i>	

Product	Item #
Ultra-Bond Plus Intro Kit	
Ultra-Bond Plus Intro Kit	032599110
<i>Contents: 1 – intro syringe shades A1, B1, B0, Supreme White, 1 – 1g syringe Try-In Paste shades A1, B1, B0, Supreme White, 3ml Etch'N'Seal, 2ml Cerinate Prime, 1g porcelain conditioner, 3g Tenure S Dab-Eze, syringe mixing tips 8-element brown/clear (10pc), syringe tips pink 18G (6pc), brush tips (25pc), brush-tip handle, mixing pad, curing rings (3pc), all-surface bonding chart, instructions</i>	
Ultra-Bond Plus Multi-Shade Kit	
Ultra-Bond Plus Multi-Shade Kit	032599750
<i>Contents: 1 – 6g UB syringe shades A1, B1, B0, Supreme White, 1 – 1g syringe Try-In Paste shades A1, B1, B0, Supreme White, 5ml Cerinate Prime, 3ml Etch'N'Seal syringe, 6g porcelain conditioner, 3g Tenure S Dab-Eze, auto-mixing tips (30pc), brush tips (25pc), syringe tips pink 18G (18pc), brush-tip handle, mixing pad, curing rings (3pc), all-surface bonding chart, instructions</i>	
Ultra-Bond Plus	
Ultra-Bond Plus Custom Kit	032599950
<i>Contents: 4 – 6g syringes, auto-mixing tips (60pc), intraoral tips (40pc), light cure test rings (3pc), instructions</i>	
Ultra-Bond Refills	
Ultra-Bond Plus Neutral Refill	032599370
Ultra-Bond Plus Supreme White Refill	032599380
Ultra-Bond Plus A1 Refill	032599310
Ultra-Bond Plus A2 Refill	032599390
Ultra-Bond Plus B0 Refill	032599350
Ultra-Bond Plus B1 Refill	032599330
<i>Contents: 1 – 6g syringe, 1 – 1g syringe Try-In Paste, auto-mixing tips (10pc), syringe tips pink 18G (6pc), curing rings (3pc), all-surface bonding chart, instructions</i>	



Nuance Universal Universal Light-Cured Restorative Material

Nuance is a light-cured micro hybrid composite that simplifies the composite restoration process like never before. It's translucent when thin and opaque when thick. One body shade can be used for an entire restoration with flawless, natural-looking results.

- Superb handling properties
- Creamy consistency easily spreads
- Available in nine shades
- Convenient flip-top syringe or pre-loaded tips for minimal waste

nuance®



Before



After



Nuance FLOW

Highly Polishable Flowable Composite

Nuance FLOW composite material is ideal for the pulpal or axial layer in posterior restorations. It has a unique self-leveling viscosity that allows for easy placement without running or dripping and it has optimal compatibility with Bis-acrylic materials.

- Easy to place
- Excellent viscosity
- Radiopaque
- Convenient syringe delivery
- Available in eight shades

Product	Item #
Nuance Universal Syringe	
Nuance Universal Syringe A1	CR1101
Nuance Universal Syringe A2	CR1102
Nuance Universal Syringe A3	CR1103
Nuance Universal Syringe A3.5	CR1104
Nuance Universal Syringe B1	CR1108
Nuance Universal Syringe B2	CR1109
Nuance Universal Syringe C3	CR1112
Nuance Universal Syringe BL (Bleach White)	CR1106
Nuance Universal Syringe OA2 (Opaque A2)	CR1121
Contents: 1 – 4g syringe, instructions	

Product	Item #
Nuance Universal Singule	
Nuance Universal Singule A1	CR1135
Nuance Universal Singule A2	CR1136
Nuance Universal Singule A3	CR1137
Nuance Universal Singule A3.5	CR1138
Nuance Universal Singule B1	CR1142
Nuance Universal Singule B2	CR1143
Nuance Universal Singule C3	CR1146
Nuance Universal Singule BL (Bleach White)	CR1140
Nuance Universal Singule OA2 (Opaque A2)	CR1155
Contents: 16 – 0.25g singules, instructions	

Product	Item #
Nuance Flowable	
Nuance FLOW A1	CR1217
Nuance FLOW A2	CR1215
Nuance FLOW A3	CR1214
Nuance FLOW A3.5	CR1213
Nuance FLOW B1	CR1212
Nuance FLOW B2	CR1211
Nuance FLOW BL1	CR1210
Nuance FLOW BL2	CR1209
Contents: 1 – 1g syringe, instructions	

Virtuoso® Universal Micro Hybrid Composite



Virtuoso Universal offers esthetics, high polishability and optimal handling all in one formula. The innovative micro hybrid technology makes Virtuoso an ideal universal restorative material for anterior and posterior applications.

- Superb polishability
- Wear and stain resistant
- Wide shade range
- Convenient pre-loaded tips for minimal waste

Product	Item #
Virtuoso Universal Kit	
Intro Kit	030381950
Contents: 10 – .25g capsules A2, A3.5 & Enamel, 3g – Gloss-N-Seal, 1.5ml Tenure Uni-Bond A, 1.5ml Tenure Uni-Bond B, 1g Virtuoso Flowable A2 syringe, micro brush, dappen dish, curing rings (3pc), syringe tips 20G (5pc), all-surface bonding chart, instructions	
Virtuoso Universal Refill Kits – 10 Capsules	
Refill Kit Shade A1	030381969
Refill Kit Shade A2	030381963
Refill Kit Shade A3	030381968
Refill Kit Shade A3.5	030381961
Refill Kit Shade B0	030381967
Refill Kit Shade B1	030381966
Refill Kit Shade B2	030381970
Refill Kit Shade C1	030381965
Refill Kit Shade C2	030381964
Refill Kit Shade – Enamel	030381962
Contents: 10 – .25g capsules, instructions	

Virtuoso® Flowable Tried-and-True Flowable Composite



Virtuoso flowable is ideal for small restorations, especially those difficult to access. Optimal handling characteristics make this light-cure, low-viscosity, flowable composite easy to use.

- Light-cure
- Low viscosity
- Precision needle tip for precise placement
- Ideal bubble-free flowable consistency
- Wide shade range

Product	Item #
Virtuoso Flowable Intro Kit	
Intro Kit	030382800
Contents: 1g Virtuoso Flowable A2, A3, B3 & C3, 3ml BondLink Bond Enhancer, syringe tip 20G (20pc), curing ring (1pc), all-surface bonding chart, instructions	
Virtuoso Flowable Refill Kits – 4 Syringes	
Refill Kit Shade A1	032381820
Refill Kit Shade A2	032381815
Refill Kit Shade A3	032381825
Refill Kit Shade A3.5	032381830
Refill Kit Shade B1	032381835
Contents: 4 – 1g syringes, instructions	
Virtuoso Flowable Custom Kit	
Flowable Custom Kit	030381800
Contents: 15 – 1g syringes of the available shades, instructions	
Virtuoso Flowable Custom Refill Kit	
Flowable Custom Refill Kit	032381810
Contents: 4 – 1g syringes of the available shades, instructions	



Before





After

Available Virtuoso Shades



Product	Item #
Ancillary Composite Items	
Singule Dispensing Gun	CR1183
SY Centrix Mark II (1pc)	38115
Syringe Tips 20G (20pc)	030007511
Syringe Tips 20G (100pc)	030007611

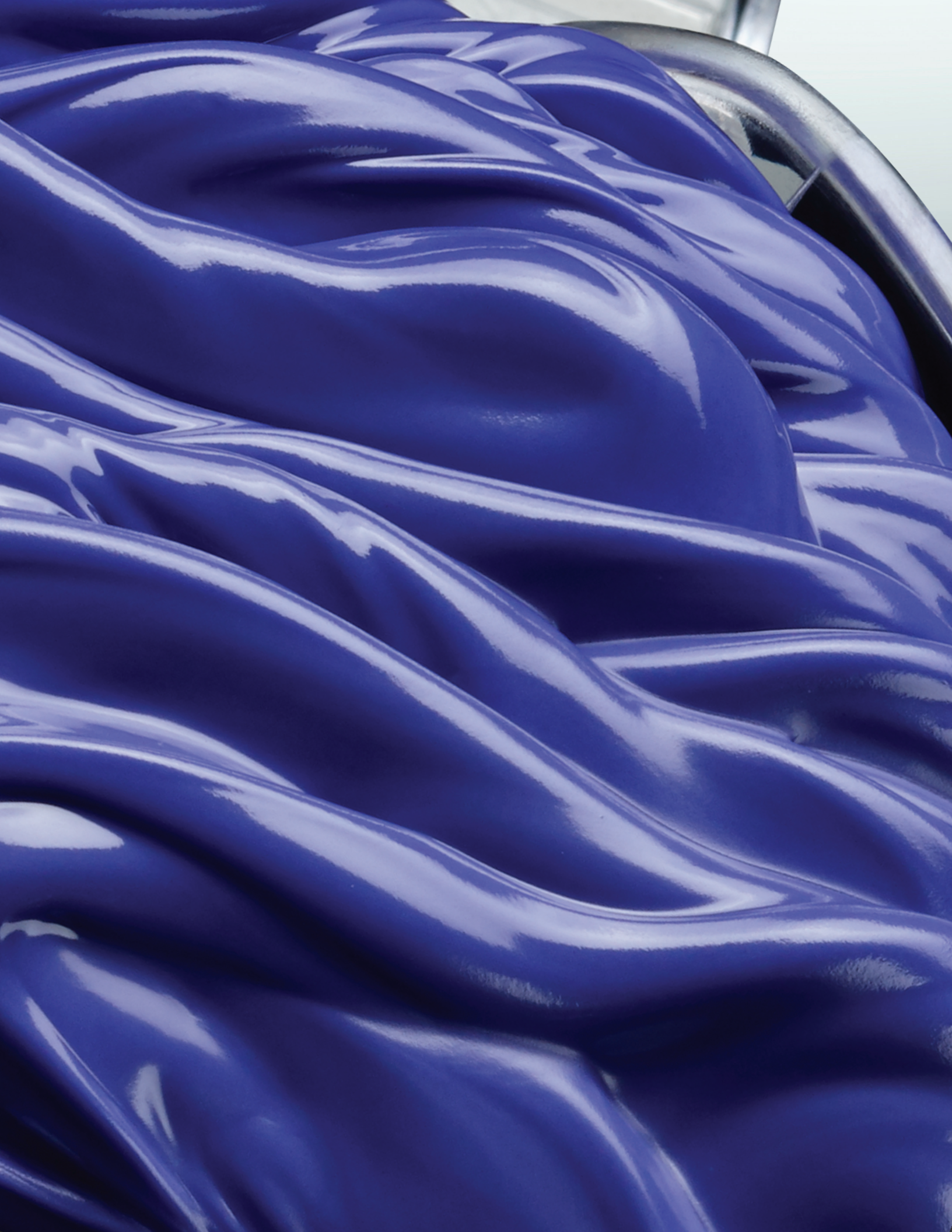

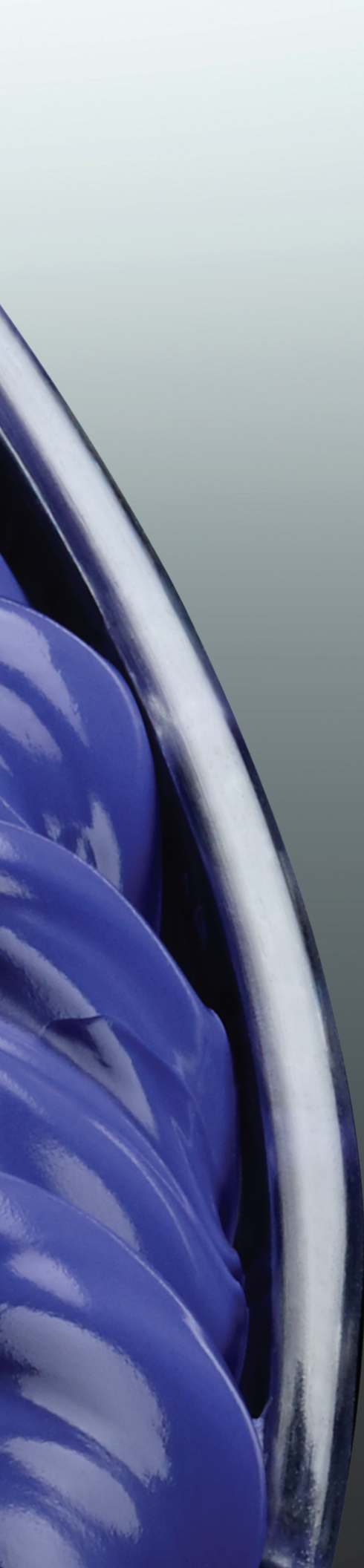

For reliable, consistent impressions, time-after-time, choose DenMat® materials. From simple cases to complex full-mouth restorations, our VPS materials and our award-winning impression trays deliver the quality you and your patients can depend on.



SplashMax

The Ultimate VPS Impression Material for Accuracy, Reliability and Value

SplashMax is a super VPS impression material that offers instant hydrophilicity, because dead-on accuracy is needed when the material first comes in contact with dentition. In independent testing, SplashMax demonstrated the lowest contact angle almost immediately.

- Extreme wettability
- Maximum tear strength and elasticity
- Excellent dimensional accuracy
- Four high-contrast viscosities
- Wild berry scent



Extra Lite



Lite



Medium



Heavy

SplashMax features the lowest contact angle within seconds!*



*Data on file.



Product	Item #
SplashMax 2-Pack Cartridges Regular Set (4:30)	
Extra Lite Body	SPD1613
Lite Body	SPD1612
Medium Body	SPD1611
Heavy Body	SPD1610
<i>Contents: 2 – 50ml cartridges, tips, instructions</i>	
SplashMax 2-Pack Cartridges Half-Time Set (2:15)	
Extra Lite Body	SPD1617
Lite Body	SPD1616
Medium Body	SPD1615
Heavy Body	SPD1614
<i>Contents: 2 – 50ml cartridges, tips, instructions</i>	
SplashMax 8-Pack Cartridges Regular Set (4:30)	
Extra Lite Body	SPD1606
Lite Body	SPD1608
<i>Contents: 8 – 50ml cartridges, instructions</i>	
SplashMax 8-Pack Cartridges Half-Time Set (2:15)	
Extra Lite Body	SPD1607
Lite Body	SPD1609
<i>Contents: 8 – 50ml cartridges, instructions</i>	
SplashMax 20-Pack Cartridges Regular Set (4:30)	
Medium Body	SPD1604
Heavy Body	SPD1602
<i>Contents: 20 – 50ml cartridges, instructions</i>	
SplashMax 20-Pack Cartridges Half-Time Set (2:15)	
Medium Body	SPD1605
Heavy Body	SPD1603
<i>Contents: 20 – 50ml cartridges, instructions</i>	
SplashMax 2-Pack Jumbo Cartridges Regular Set (4:30)	
Medium Body	SPD1601
Heavy Body	SPD1600
<i>Contents: 2 – 360ml jumbo cartridges, mixing tips, rotating rings, instructions</i>	

To order, call 1-800-4DENMAT (1-800-433-6628) or visit www.denmat.com

dm
DenMat®




Splash!

A Hydrophilic VPS Impression Material You Can Count On

Splash! is the ideal choice for clinicians who demand reliability and value in a VPS impression material. It features excellent flow properties, tear strength and a scent that patients love.

- Super hydrophilic
- Two set times (2:15 and 4:30)
- Available in five high-contrast colors and four delivery methods
- Wild berry scent

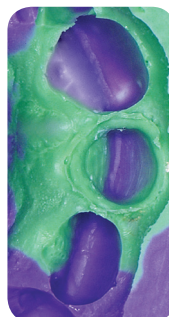

5 great tips for taking accurate impressions with *all* our impression materials!

- 1 Line plastic trays with an adhesive before extruding the material into the tray to prevent separation of the impression from the tray.
- 2 Avoid overfilling trays with impression material — only fill to the top of the tray and limit the amount of material extruded onto the pallet (for uppers). This saves material and prevents patient from gagging.
- 3 Use intraoral tips when applying the light body wash to ensure the material flows into all areas of the preparation.
- 4 Seat the impression on the front and back teeth simultaneously to ensure ultimate accuracy.
- 5 Let the impression set for the full duration of working and setting time to avoid the risk of pulls or tears. Using a timer to prevent removal of the impression too early is always a good idea.

Product	Item #
Splash! 2-Pack Cartridges Regular Set (4:30)	
Extra Lite Body	SPD1208
Lite Body	SPD1209
Medium Body	SPD1210
Heavy Body	SPD1211
<i>Contents: 2 – 48ml cartridges, tips, instructions</i>	
Splash! 2-Pack Cartridges Half-Time Set (2:15)	
Extra Lite Body	SPD1230
Lite Body	SPD1231
Medium Body	SPD1232
Heavy Body	SPD1233
<i>Contents: 2 – 48ml cartridges, tips, instructions</i>	
Splash! Bulk Packs Regular Set (4:30)	
Extra Lite Body (8-pack)	SPD1587
Lite Body (8-pack)	SPD1588
Medium Body (20-pack)	SPD1589
Heavy Body (20-pack)	SPD1590
<i>Contents: Extra Lite & Lite: 8 – 48ml cartridges; Medium & Heavy: 20 – 48ml cartridges</i>	
Splash! Bulk Packs Half-Time Set (2:15)	
Extra Lite Body (8-pack)	SPD1591
Lite Body (8-pack)	SPD1592
Medium Body (20-pack)	SPD1593
Heavy Body (20-pack)	SPD1594
<i>Contents: Extra Lite & Lite: 8 – 48ml cartridges; Medium & Heavy: 20 – 48ml cartridges</i>	
Splash! FatPak Intro Kits Regular Set (4:30)	
Medium Body	SPD1532
Heavy Body	SPD1531
<i>Contents: 1 – 300ml base, 1 – 62ml catalyst, cartridge body, tips, instructions</i>	
Splash! FatPak Refill Kits Regular Set (4:30)	
Medium Body	SPD1530
Heavy Body	SPD1529
<i>Contents: 2 – 300ml base, 2 – 62ml catalyst, tips, instructions</i>	
Splash! Putty Jars	
Regular Set Putty Jars (5:30)	SPD1212
Half-Time Set Putty Jars (2:45)	SPD1235
<i>Contents: 1 – 250ml jar base, 1 – 250ml jar catalyst, mixing spoons, instructions</i>	
Splash! Putty Paks	
Regular Set Putty Paks (5:30)	SPD1219
Half-Time Set Putty Paks (2:45)	SPD1237
<i>Contents: 12 – 18ml paks, instructions</i>	

Product	Item #
Ancillary Impression Items	
Mixing Tips Yellow (50pc)	SPD1182
Mixing Tips Blue (50pc)	SPD1183
Intraoral Tips Yellow (50pc)	SPD1185
FatPak Dynamic Mixers – <i>Splash & Precision Only</i> (50pc)	SPD1528
Jumbo Cartridge Mixing Tips & Rotating Ring – <i>SplashMax Only</i> (11pc)	030897101
Cartridge Dispenser Gun 50ml – 1:1	030799105



Precision


Exceptional Detail Reproduction and Moisture Displacement the First Time... and Every Time

Precision is a premium VPS impression material with amazing detail reproduction that displaces moisture to create highly detailed margins in even the most challenging situations.


- Outstanding dimensional stability
- Thermally activated flexible set times
- Wild berry scent
- Available in five high-contrast colors and four delivery methods



Product	Item #
Precision 2-Pack Cartridges Flex Set Time (2:45 – 3:45)	
Extra Lite Body	SPD1508
Lite Body	SPD1509
Medium Body	SPD1510
Heavy Body	SPD1511
Contents: 2 – 48ml cartridges, tips, instructions	
Precision FatPak Intro Kits Flex Set Time (2:45 – 3:45)	
Medium Body	SPD1536
Heavy Body	SPD1535
Contents: 1 – 300ml base, 1 – 62ml catalyst, cartridge body, tips, instructions	
Precision FatPak Refill Kits Flex Set Time (2:45 – 3:45)	
Medium Body	SPD1534
Heavy Body	SPD1533
Contents: 2 – 300ml base, 2 – 62ml catalyst, tips, instructions	
Precision Putty Jars Flex Set Time (2:45 – 3:45)	
Putty Jars	SPD1517
Contents: 1 – 250ml jar base, 1 – 250ml jar catalyst, mixing spoons, instructions	
Precision Putty Paks Flex Set Time (2:45 – 3:45)	
Putty Paks	SPD1514
Contents: 12 – 18ml paks, instructions	

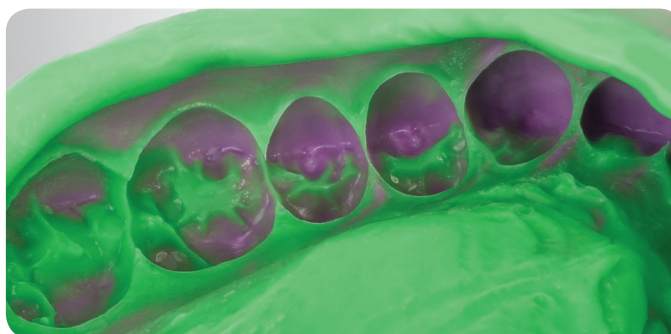


Traditional Impression Material

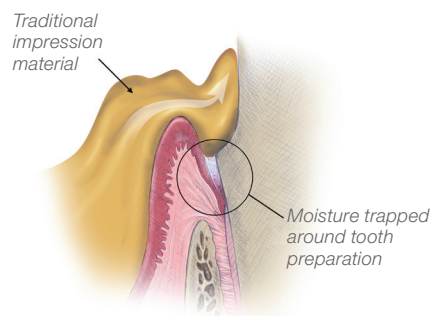

Traditional VPS material fails to achieve complete moisture displacement

Precision VPS Impression Material

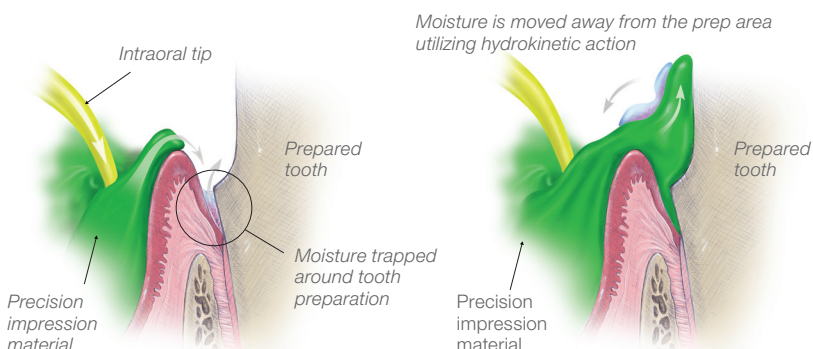

Hydrokinetics in Action: Seeking Out & Repelling Moisture
(Detail of Precision material entering tooth sulcus)



Product	Item #
Bite Registration	
Vanilla Bite	SPD1105
Chocolate Bite	SPD1186
Clear Bite	SPD1104
MegaBite	SPD1370
<i>Contents: 2 – 48ml cartridges, tips, instructions</i>	



Bite Registration Materials

Stability, Accuracy and Rigidity... Everything You Need for Bite Registrations

DenMat offers a complete line of VPS bite registration materials to meet all of your procedural needs. Easy to extrude and trim, these materials are used for closed, open- and double-bite registrations.

- Extreme accuracy and stability
- Fast set times
- Easy trimming

Mousse-like Consistency



Vanilla Bite™

- Rigid
- Quick set (0:55)
- Ideal for 1 – 2 unit cases
- Extremely accurate and dimensionally stable

Medium Rigidity



Chocolate Bite™

- Slower set time (1:30)
- Ideal for taking full-arch quadrant registrations
- Allows for longer working times

Transparent & Flexible



Clear Bite™

- Transparent
- Ideal for temporary matrix (2:30 set time)
- Enables light passage for curing
- Great visibility of voids or bubbles

Super Rigid



MegaBite®

- Very rigid VPS bite registration material
- Fast set time (0:55)
- Exceptional accuracy
- Ideal for centric occlusion and centric relation registrations
- Minimized occlusal adjustments



Product	Item #
The Gripper	
Anterior	SPD1414
Posterior	SPD1415
3/4-Arch	SPD1416
Full-Arch Large	SPD1417
Full-Arch	SPD1418
Sideless Posterior	SPD1541
<i>Contents: 40 trays per box</i>	



Anterior



Posterior



3/4-Arch



Full-Arch




Sideless Posterior

The Gripper® Super Rigid Impression Trays

The super rigid design of The Gripper makes it ideal for all VPS impression applications.

- Dual-arch disposable tray
- Rigid side walls and soft-lined mesh for exceptional accuracy
- Ergonomic shape for superior placement and comfort

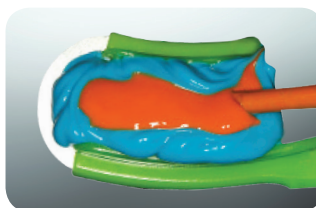

Full-Arch

Anterior

Quad-Arch

Perfectray®


Our Award-Winning Line of Trays Offer the Perfect Tray for Any Procedure

Perfectrays are crystal-clear acrylic impression trays available in many sizes for every patient and every impression. They are super rigid with ergonomic handles, finger rests and perforated for superior material retention. With nine different options in a variety of sizes and shapes, you're sure to find the perfect tray for any procedure.

- Single-arch disposable tray
- Perforations provide superior retention
- Available in a variety of sizes

Product	Item #
Perfectray	
#1 Large Upper	SPD1400
#2 Large Lower	SPD1401
#3 Medium Upper	SPD1402
#4 Medium Lower	SPD1403
#5 Small Upper	SPD1404
#6 Small Lower	SPD1405
#7 Quad. Upper R/Lower L	SPD1406
#8 Quad. Upper L/Lower R	SPD1407
#9 Anterior	SPD1408
<i>Contents: 20 trays per box</i>	



A close-up photograph of a person's mouth, showing a bright smile with white teeth. The image is partially obscured by a dark grey, angular graphic element in the lower right corner, which contains white text. The background of the image is a soft, out-of-focus gradient of light blue and white.

DenMat® has a long history of innovation in whitening. From the launch of the original Rembrandt® Whitening in 1986 through Lumibrite® and LumiSmile® White, DenMat has always been at the forefront of safe, effective teeth whitening.



LumiSmile White
Take-Home Bulk Pack



LumiSmile White
Take-Home Starter Kit



LumiSmile White
Take-Home Refill Kit



Grow your
practice with
LumiSmile White's
high-impact
marketing materials.
[View online now!](#)

LumiSmile White Take-Home Professional Take-Home Whitening System

LumiSmile White's proprietary formula offers premium professional-strength whitening for every smile.* LumiSmile White will whiten your patient's smile in 1 – 2 weeks with minimal sensitivity. Choose either 32%, 22% or 16% carbamide peroxide concentrations flavored with real peppermint oil. LumiSmile White features an 18-month shelf life with no refrigeration required.† Brighten your practice today with LumiSmile White and get more patients to say YES to a whiter, brighter, more confident smile.

Product	Item #
LumiSmile White Take-Home Kits†	
16% Mint	WHT1016
22% Mint	WHT1022
32% Mint	WHT1032
16% Mint No-Fluoride	WHT1016NF
22% Mint No-Fluoride	WHT1022NF
32% Mint No-Fluoride	WHT1032NF
<i>Contents: 6 syringes (18 applications), 2 resin sheets, mouthguard case, instructions</i>	
LumiSmile White Take-Home Refill Kits†	
16% Mint	WHT1116
22% Mint	WHT1122
32% Mint	WHT1132
16% Mint No-Fluoride	WHT1116NF
22% Mint No-Fluoride	WHT1122NF
32% Mint No-Fluoride	WHT1132NF
<i>Contents: 12 syringes (36 applications)</i>	
LumiSmile White Take-Home Bulk Pack†	
22% Mint Non-Personalized	WHT5022
32% Mint Non-Personalized	WHT5032
22% Mint Personalized	WHT5022P
32% Mint Personalized	WHT5032P
<i>50 individually wrapped syringes with patient instructions, 1 practice building guide, 5 interchangeable marketing postcards, 25 promotional cards</i>	

*Individual results may vary.

†Store at or below 75° F. Do not freeze.

‡Some products may not be available in Canada.



LumiSmile White In-Office*

Professional In-Office Whitening System

LumiSmile White In-Office whitening system is a premium chairside whitening solution for every smile*. Formulated to provide effective and fast whitening with minimal sensitivity, its 25% hydrogen peroxide concentration gives patients an effective and enjoyable whitening experience in just an hour. For optimal results, patients can maintain their bright white smile at home with LumiSmile White Take-Home.

- 25% hydrogen peroxide
- Minimal sensitivity
- Long unrefrigerated shelf life

Product	Item #
LumiSmile White In-Office Kits*	
LumiSmile White In-Office Kit	WHT0025
<i>Contents: 1 whitening syringe (1 application), 1 desensitizing gel syringe, 1 paint-on dental dam syringe, lip moisturizer, mini brush applicator, syringe tip 19G, lip expander, dappen dishes (3pc), instructions</i>	
LumiSmile White In-Office & Take-Home Combo Kit	WHT0030
<i>Contents: 1 in-office whitening syringe (1 application), 1 desensitizing gel syringe, 1 paint-on dental dam syringe, lip moisturizer, mini brush applicator, syringe tip 19G, lip expander, dappen dishes (3pc), 2 - 2.5ml take-home syringes (22% carbamide peroxide), instructions</i>	



Lumibrite® Take-Home

Professional Take-Home Teeth Whitening

Lumibrite Take-Home is gentle, efficacious and proven to whiten smiles in 1 – 2 weeks. Choose either 32%, 22% or 16% carbamide peroxide concentrations in mint or fruit flavors to meet the specific needs and preferences of your patients.

- 32%, 22% and 16% carbamide peroxide
- Choice of mint or tropical fruit flavors
- Available in:
 - Standard Kits (6 syringes, 18 applications)
 - Refill Kits (12 syringes, 36 applications)

Product	Item #
Lumibrite Take-Home Kits	
16% Mint	034359145
22% Mint	034339145
32% Mint	034139145
16% Fruit	034369145
22% Fruit	034349145
32% Fruit	034239145
<i>Contents: 6 syringes (18 applications), 2 resin sheets, mouthguard case, instructions</i>	
Lumibrite Take-Home Refill Kits	
16% Mint	034358725
22% Mint	034338725
32% Mint	034138725
16% Fruit	034368725
22% Fruit	034348725
32% Fruit	034238725
<i>Contents: 12 syringes (36 applications)</i>	




Lumibrite Chairside*


Professional Chairside Teeth Whitening

The Lumibrite Chairside Whitening System is ideal for fast, effective and gentle whitening with minimal sensitivity. Lumibrite Chairside Whitening is effective with or without a light.

- 35% hydrogen peroxide
- Includes paint-on dental dam and accessories
- Unbeatable value compared to the leading competition

Product	Item #
Lumibrite Chairside Whitening System*	
Chairside Whitening System (US)	035359155
Chairside Whitening System (Canada)	035359159
<i>Contents: 1 whitening syringe (1 application), 1 desensitizing gel syringe, 1 paint-on dental dam syringe, lip moisturizer, mini brush applicator, syringe tip 19G, lip expander, dappen dishes (3pc), instructions</i>	
Lumibrite Desensitizing Kit	
Desensitizing Kit	034159555
<i>Contents: Desensitizing gel syringe (3 applications)</i>	

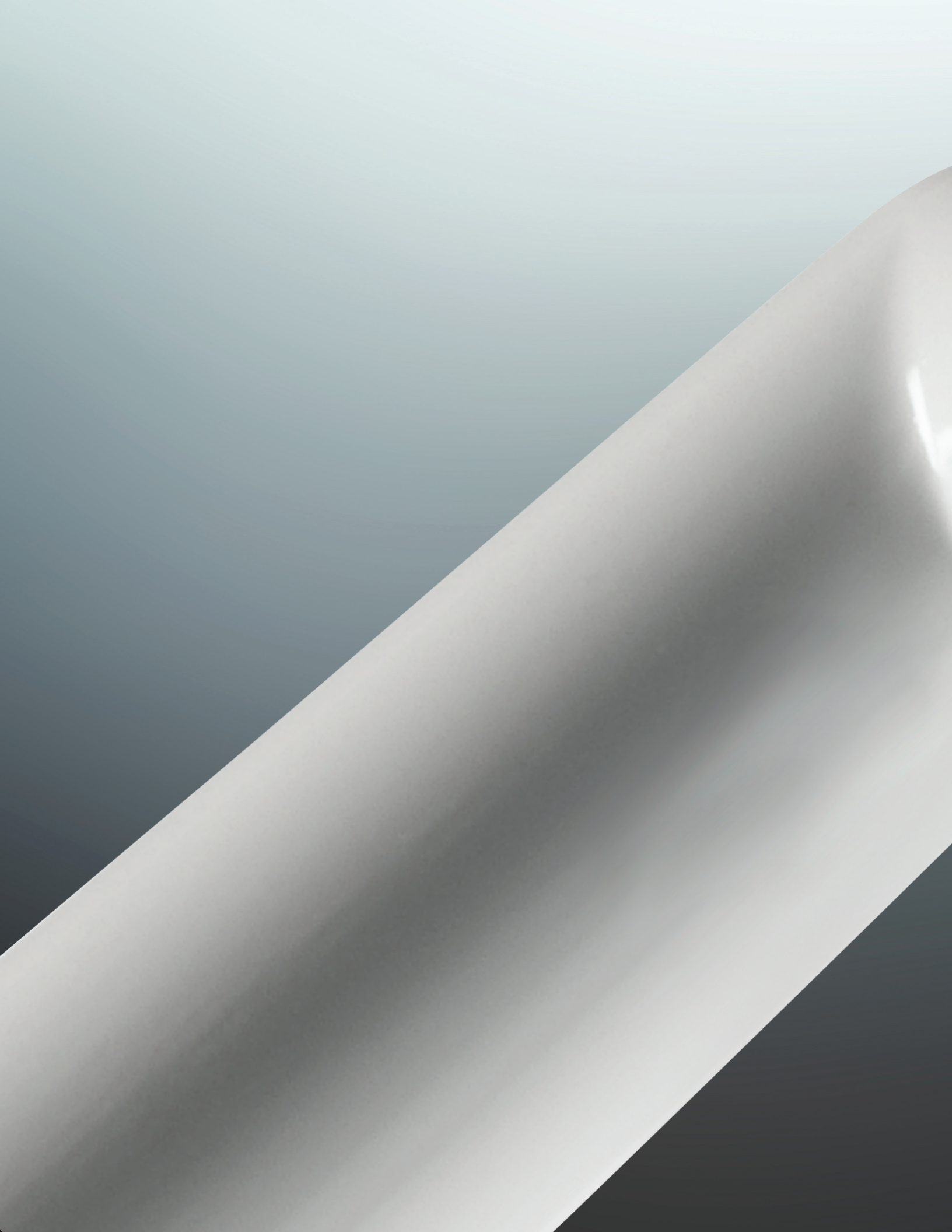

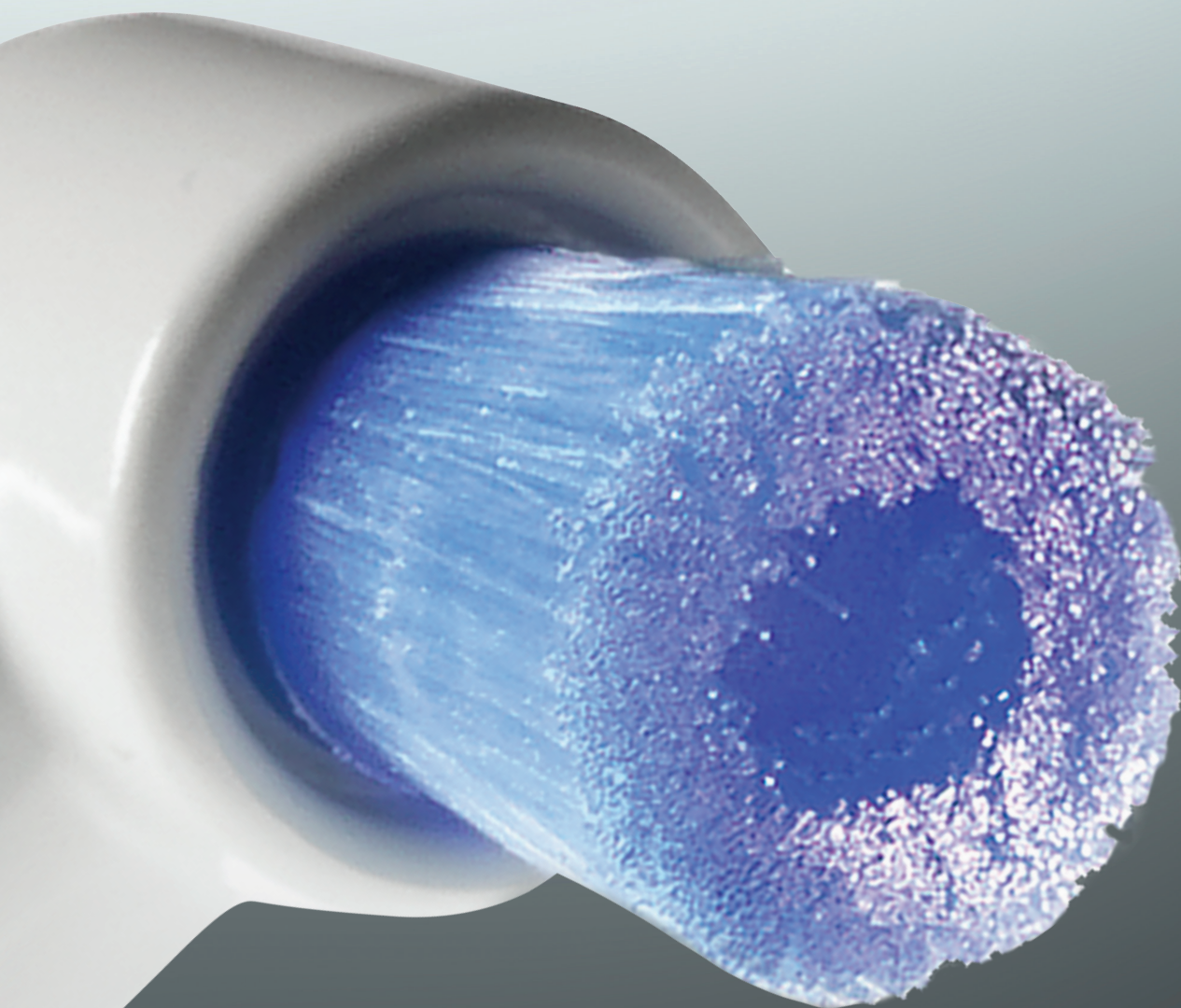

DenMat® offers a complete array of oral hygiene products to meet the needs of you and your patients. From our ultrasonic scaling system to our power toothbrush that allows patients to maintain their regimen at home, we've got you covered.

Whether in your chair or at home, offer your patients superior clinical care.

POWER TOOTHBRUSHES



Rotadent ProCare

The Professional Rotary Toothbrush

For over 25 years, the Rotadent power toothbrush has established itself as the premiere at-home solution for periodontal maintenance. DenMat is proud to introduce the latest generation of the Rotadent toothbrush: Rotadent ProCare.

Rotadent ProCare is a true extension of in-office procedures, helping you ensure that your patients maintain and improve their oral health between office visits. It's also been shown to be a "one-step" system, which simplifies and streamlines homecare for patients. Here are just a few of the proven facts from clinical studies:

Rotadent...

- is the only power toothbrush clinically proven to be as effective as brushing and flossing combined¹
- removes 92% of plaque in the first minute of brushing²
- consistently out-performed other brushes in in-vitro stain removal from both flat and interproximal surfaces³
- increases patient compliance with homecare and recare appointments⁴

Unlike competitive power toothbrushes, Rotadent ProCare is available to patients **EXCLUSIVELY** via dental professionals.


Rotadent
ProCare™
The Professional Rotary Toothbrush



How does it work?

Proprietary MicroFilaments™

Each Rotadent brush head has **4,577** proprietary MicroFilaments that mechanically sweep away plaque and biofilm. Each MicroFilament is 1/3 the width of a traditional toothbrush bristle.

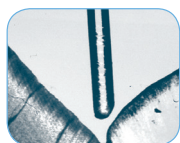


Ordinary Bristles



Rotadent MicroFilaments



Rotadent ProCare MicroFilaments effectively reach below the gumline and interproximally.

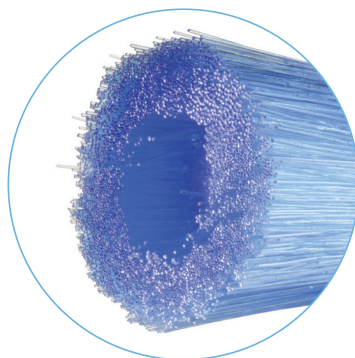

Unique Handle Design

The Rotadent toothbrush was designed to be an at-home extension of the Prophy handpiece that is used to fight plaque and control disease during Hygiene appointments.



3 Brush Head Options



Hollow-Tip

- Designed for cleaning and polishing surfaces of the teeth and sweeping under the gumline
- Flat, cup-like design makes it ideal for absorbing and applying oral medication



Short-Tip

- Effective for cleaning and polishing surfaces of the teeth
- Able to deliver oral medication more deeply into periodontal pockets and between teeth



Long-Tip

- Most effective for targeting hard-to-reach areas
- Sweeps plaque and debris away from braces, multi-unit fixed bridges or between teeth

Who do you dispense to?

Rotadent is not just a solution for periodontal patients. It is truly optimized to treat all patient types.

Orthodontic

- Controls gingivitis and plaque that build up around braces⁵
- Helps minimize staining and decalcification that many Ortho patients experience⁶

Implant

- Safe home oral hygiene regimen to use on implant surfaces⁷
- Reduces surface abrasion, which may lead to reduced buildup of plaque⁷

Esthetic

- Uses less brushing force than competitive brushes⁸
- Ideal solution for veneer patients with extra-sensitive soft tissue

Restorative

- Safe for use on dental restorations, including crowns, fillings and bridges
- Effectively balances low abrasivity with high-cleaning efficiency³



Product	Item #
Rotadent ProCare*	
Rotadent ProCare	94000N
Rotadent ProCare Staff Units	94000S
<i>Contents: power toothbrush, charger, (2) brush heads: (1) hollow, (1) short</i>	
Rotadent ProCare/Contour Brush Heads	
Hollow-Tip (1-pack)	9290H
Short-Tip (1-pack)	9290S
Long-Tip (1-pack)	9290L
Rotadent Plus Brush Heads	
Hollow-Tip (2-pack)	9094H
Short-Tip (2-pack)	9094S
Long-Tip (2-pack)	9094L
Rotadent Classic Brush Heads	
Black Hollow-Tip (1-pack)	9090H
Black Short-Tip (1-pack)	9090S
Black Long-Tip (1-pack)	9090L

*2-year warranty

1. Boyd RL, Murray P, Robertson PB. Effect on periodontal status of rotary electric toothbrushes vs. manual toothbrushes during periodontal maintenance. *J Periodontol.* 1989;60(7):390-395.
 2. Preber H, Ylipää V, Bergström J, Rydén H. A comparative study of plaque removing efficiency using rotary electric and manual toothbrushes. *Swed Dent J.* 1991;15(5):229-234.
 3. Sarker S, McLay L, Boyd RL. Clinical and laboratory evaluation of powered electric toothbrushes: laboratory determination of relative interproximal cleaning efficiency of four powered toothbrushes. *Clin Dent.* 1997;8(3):81-85.
 4. Bader HL. Compliance with a rotary powered brush: a 15-year review of patient compliance. *RDH.* May 2011; 100 - 105.
 5. Boyd RL, Murray P, Robertson PB. Effect of rotary electric toothbrush versus manual toothbrush on periodontal status during orthodontic treatment. *Am J Orthod Dentofacial Orthop.* 1989;96(4):342-347.
 6. Boyd RL, Rose CM. Effect of rotary electric toothbrush versus manual toothbrush on decalcification during orthodontic treatment. *Am J Orthod Dentofacial Orthop.* 1994;105(5):450-456.
 7. Thomson-Neal D, Evans GH, Meffert RM. Effects of various prophylactic treatments on titanium, sapphire, and hydroxyapatite-coated implants: an SEM study. *Int J Periodontics Dent.* 1989;9(4):300-311.
 8. Boyd RL, McLay L, Zahradnik R. Clinical and laboratory evaluation of powered electric toothbrushes: in vivo determination of average force for use of manual and powered toothbrushes. *J Clin Dent.* 1997; 8(3 Spec. No): 72 - 75.

To order, call **1-800-4DENMAT** (1-800-433-6628) or visit www.denmat.com

dm
DenMat®

ULTRASONIC SCALERS



PSP® Pro-Select
Platinum®
PIEZO ULTRASONIC

PSP Piezo Ultrasonic Scaler

Two instruments in one — ultrasonic scaling and heated irrigation system. PSP is a superb piezo ultrasonic scaling solution and the only system that combines scaling with heated subgingival irrigation. Upgrade your practice today.

PSP — Piezo Ultrasonic Scaler

Easy on you. Easy on your patients.

- Patient comfort.**
 Patient acceptance is increased because of PSP's high-frequency linear tip motion, quiet operation, minimal water output, and effectiveness with light, lateral pressure.
- Clinician comfort.**
 The lightweight, ergonomic handpiece is designed to reduce muscle strain. Less pressure is required to effectively scale. The clinician can work more efficiently and can perform gross debridement and definitive root planing in one step, without changing tips.
- Closed system.**
 Multi-fluid system offers choice of using distilled water, chlorhexidine and fluoride in a simple, three medicament bottle design. Because the fluid delivery system is closed, infection control is enhanced and the risk of biofilm contamination is reduced.



Product	Item #
PSP Standard Package	
PSP Standard Package	S42350
<i>Contents: 1 – PSP control unit, 1 – scaler handpiece, 2 – LED handpiece lights, 1 – irrigation handpiece, 2 – autoclavable cassettes, 3 – universal straight tips, 2 – six-foot autoclavable cables, 2 – tip change kits, 10 – cannulas, 6 – solution bottles, 4 – 16oz bottles of chlorhexidine, 1 – 64oz bottle of 2.0% NaF irrigation solution* 1 – multi-function foot control, 1 – power transformer, 1 – operating manual, 1 – one-year warranty</i>	
PSP High-Production Package	
PSP High-Production Package	S42352
<i>Contents: All items in the standard package plus 2 additional scaler handpieces, 1 additional LED handpiece light and 1 additional 6-foot autoclavable cable</i>	
PSP Accessories	
Scaler Handpiece	S40801
Irrigation Handpiece	S40851
Universal Straight Tip	S30930
Universal Left-Bend Tip	S30940
Universal Right-Bend Tip	S30950
Set of 3 Straight Tips	S30960
Set of 3 – 1 Straight, 1 Right-Bend, 1 Left-Bend	S30920
Tip Grip	100540
Tip Gauge/Wrench	100773
PSP Replacement Parts	
LED Light Attachment	S40871
Irrigation Cannulas (100/box)	2054
Transformer/Adapter	S30500
Replacement Irrigation Bottle (1)	S40200
6-ft Autoclavable LED Cable	S40750
Foot Pedal w/ Cable	S40901

Gum Disease Symptoms

Healthy



Signs of a healthy mouth:

- Pink gum color
- Gumline hugs teeth tightly
- No bleeding

Gingivitis



Symptoms of gingivitis:

- Gums bleed
- Inflamed red gum tissue
- Bad breath and bad taste

Early Periodontitis



Symptoms of early periodontitis:

- Gums may pull away from teeth
- Inflamed red gum tissue
- Bad breath and bad taste
- Bone loss around teeth
- Pockets of bacterial infection of 4mm between teeth and gums

Moderate Periodontitis



Symptoms of moderate periodontitis:

- Periodontal abscesses
- Gums recede and teeth look longer
- Loose or separating teeth
- Bad breath and bad taste
- Pockets of bacterial infection from 4 to 6mm deep

Advanced Periodontitis



Symptoms of advanced periodontitis:

- Mobile or loose teeth
- Persistent bad breath and bad taste
- Exposed roots and sensitive to hot and cold
- Loss of the adjacent supporting bone
- Pockets of bacterial infection in excess of 6mm deep

PREVENTATIVES



Pro-DenRx 2.0% Neutral Sodium Fluoride Rinse* In-Office Fluoride Rinse

Comes in easy-to-use single bottle with pre-measured pump to eliminate messy mixing. Safe to use on patients with restorations. Available in 64oz bottles. Mint or berry fresh flavor.

Pro-DenRx 0.2% Neutral Sodium Fluoride Rinse* Homecare Fluoride Rinse

Homecare fluoride rinse for caries prevention with plaque-controlling formula. Ideal for cosmetic patients. Available in 16oz bottles. Mint flavor.

Pro-DenRx 0.63% Stannous Fluoride Rinse* Homecare Fluoride Rinse

Homecare fluoride rinse that may prevent decay, inhibit plaque and reduce sensitivity in gums. Can be used as a rinse or with an irrigator. Available in 10oz bottles. Mint flavor.

Pro-DenRx 0.044% Daily APF Rinse* Homecare Fluoride Rinse

Homecare rinse with high fluoride uptake. Great option for patients undergoing orthodontic treatments. Available in 16oz bottles. Mint flavor.



*Pro-DenRx
2.0% Neutral Sodium
Fluoride Rinse
(In-Office)*



*Pro-DenRx
0.2% Neutral Sodium
Fluoride Rinse
(Homecare)*



*Pro-DenRx
0.63% Stannous
Fluoride Rinse
(Homecare)*



*Pro-DenRx
.044% Daily
APF Rinse
(Homecare)*



Pro-DenRx Fluoride Varnish* 5% Sodium Fluoride Varnish

Pro-DenRx varnishes come in single unit doses eliminating the need for clean-up or mixing. All three flavors are sweetened with Xylitol, a natural caries inhibitor.[†] Application provides 22,600 ppm of fluoride that may offer hypersensitivity relief for up to six months. Varnish sets rapidly in the presence of saliva, resulting in improved patient comfort and acceptance. Bubblegum, mint or melon flavor.



Pro-DenRx Fluoride Gels & Pastes*

1.1% Neutral Sodium Fluoride Gel

Prescription strength gel offers high-potency topical fluoride for caries and plaque control. Available in 2oz tubes. Mint flavor.

1.1% Plus Neutral Sodium Fluoride Dentifrice Paste

Prescription strength dentifrice combines high-potency fluoride with low-abrasive silica for a one-step brushing treatment. Available in 2oz tubes. Mint or cherry limeade flavor.



Pro-DenRx Chlorhexidine Gluconate* Antimicrobial Oral Rinse

0.12% chlorhexidine gluconate oral rinse offers antimicrobial benefits for use between dental visits for the treatment of gingivitis. Recommended for use in office or as part of a home treatment program. Available in 16oz bottles. Mint flavor.

Product	Item #
Pro-DenRx 2.0% Neutral Sodium Fluoride Rinse*	
Mint Flavor	2037MTDT
Berry Fresh Flavor	2037RBDT
Pro-DenRx 0.2% Neutral Sodium Fluoride Rinse*	
Mint Flavor	2107MTM
Pro-DenRx 0.63% Stannous Fluoride Rinse*	
Mint Flavor	2106MTDT
Pro-DenRx .044% Daily APF Rinse*	
Mint Flavor	2108MTD

Product	Item #
Pro-DenRx Fluoride Varnish Trial Pack*	
Bubblegum Flavor	VN035BG
Mint Flavor	VN035MT
Melon Flavor	VN035ML
<i>Contents: 35 – 0.5ml unit doses</i>	
Pro-DenRx Fluoride Varnish Bulk Pack*	
Bubblegum Flavor	VN200BG
Mint Flavor	VN200MT
Melon Flavor	VN200ML
<i>Contents: 200 – 0.5ml unit doses</i>	

Product	Item #
Pro-DenRx 1.1% Neutral Sodium Fluoride Gel*	
Mint Flavor	2240MTM
Pro-DenRx 1.1% Neutral Sodium Fluoride Dentifrice Paste*	
Mint Flavor	2250MTM
Cherry Limeade Flavor	2250CLM

Product	Item #
Pro-DenRx Chlorhexidine Gluconate*	
Mint Flavor	2193B



Why ATRIDOX?

- ✓ Reduces probing depth
- ✓ Improves clinical attachment
- ✓ Less bleeding upon probing

ATRIDOX®
(doxycycline hyclate) 10%

ATRIDOX

Focused Periodontal Disease Treatment

ATRIDOX (doxycycline hyclate) 10% is a major advancement in the treatment of periodontal disease. It is an easy-to-use, cost-effective, locally applied antibiotic (LAA) that is placed gently below the gum line into periodontal pockets where bacteria thrive.

ATRIDOX is the *only* LAA that is clinically proven to achieve all three of the following outcomes in treating adult chronic periodontitis:

- Gain in clinical attachment*
- Reduction of probing depths*
- Reduction of bleeding on probing*

ATRIDOX Key Features

- 1 Cost per site treated is less expensive than leading competitors**
- 2 Flows freely and easily to the bottom of pocket to adapt to root morphology
- 3 Controlled release of doxycycline for a period of 21 days*
- 4 Bioabsorbable — no removal required
- 5 Allows precise placement for targeted therapy
- 6 Dosing flexibility allows for treatment of multiple sites with a single syringe



ATRIDOX is the only LAA gel that flows deep into infected pockets

Product	Item #
ATRIDOX	
ATRIDOX (U.S.)	T10629
ATRIDOX (Canada)	T10656
Contents: 6 syringes	

*Data on file. **After 3 sites.



ViziLite® TBlue® Annual Oral Cancer Screening

ViziLite TBlue is designed to accurately detect abnormal lesions and rule out trauma and false positives. You are the only person who can spot this disease early enough to give your patients a chance at survival and a high quality of life after diagnosis. Have you helped save a life today?

- Able to reach all areas of the mouth
- Helps prevent false positives
- Improves chances of diagnosing oral cancer early

Product	Item #
ViziLite TBlue	
ViziLite TBlue (10-pack)	650199
Eyewear	650101
TBlue Supplemental Pack (6 marking systems)	650106
Angled Exam Mirrors (6-pack)	650119



RotaPoint® Interdental Cleaning Device

Made of unbreakable polypropylene, these cleaners are excellent for removing dental plaque and debris from interproximal spaces. Highly flexible, RotaPoints can easily be bent or shaped before or after insertion between teeth. RotaPoints are available in white and dispensers can be personalized with your dental practice information.

Product	Item #
Pro-DenRx RotaPoint	
RotaPoint (100, 30-count boxes)	NPRPZ
RotaPoint - Personalized* (100, 30-count boxes)	PRPZ

*All personalized items include office name, contact information and logo. A \$10 setup fee is required per order.

Perio probes have high-contrast marking for clear visibility



Sensor probes alert clinicians when universal pressure has been applied

Pro-DenRx Probes

Perio Probes

Designed for increased sensitivity, accuracy and patient-friendly comfort, these probes incorporate the latest thermoplastic technology. The black markings provide high-contrast for more accurate pocket measurements, especially in the distal regions of the oral cavity.

Sensor Probes

These innovative probes make it possible to maintain consistency of probing force between clinicians and examinations.

Product	Item #
Pro-DenRx Perio Probes (12-pack)	
CPITN Type: — 2 • 3 • 4 • 5 • 7 • 9	7101
Type: CP-12 — 3 • 6 • 9 • 12	7107
Type: 2M-12 — 2 • 4 • 6 • 8 • 10 • 12	7111
ProDenRx Sensor Probes (12-pack)	
CPITN Type: — 2 • 3 • 4 • 5 • 7 • 9	7113
Type: 2M-12 — 2 • 4 • 6 • 8 • 10 • 12	7123
CPITN Type: C — 3.5 • 5.5 • 8.5 • 11.5	7115
Type: CP-12 — 3 • 6 • 9 • 12	7119



Pro-Ties® Instrument Bundling System

Get organized with Pro-Ties, an efficient bundling system for instrument sterilization. Pro-Ties are autoclavable, stretch up to 300% and are available in four colors.


Product	Item #
Pro-DenRx Pro-Ties (6-pack)	
Yellow	650102Y
Green	650102G
Red	650102R
Blue	650102B





DenMat® proudly offers Hartzell® dental and surgical instruments. Carefully fashioned by skilled craftspersons, our instruments have maintained the same high standards of manufacturing excellence established by their founders in 1935.

Only the highest-quality components and production methods are used in manufacturing our instruments, yielding precision and consistency unparalleled in the industry.



HYGIENE INSTRUMENTS



Hygiene Instruments

Precision Hygiene Instruments

All of Hartzell's hygiene instruments are manufactured, polished and finished in the USA. With consistent bends and second-to-none edge-retention, these precision instruments will become a standard part of your perio armamentarium. Use them during routine prophylaxis appointments, or for scaling and root planing.



*Instruments are delivered pre-filled in the cassette and are not replaceable. Above photo is not representative of the assortment of instruments in the HZ-SRPKIT. Additional instruments to round out your armamentarium can be found on the opposing page.

Product	Item #
Scaling and Root Planing Kit*	
Scaling and Root Planing Kit*	HZ-SRPKIT
Contents*	
Cassette (10-Instrument)	IC10
Mirror	M5EA
Mirror Handle	CSH2A
Explorer	XA11/12
Probe (UNC 12mm)	PUNC12
Sickle Scaler	SAMJ1/2NSS
Universal Curette	SA4L/4RSS
Rigid Gracey Curette #5/6	SAGL5/6SS
Rigid Gracey Curette #7/8	SAGL7/8SS
Rigid Gracey Curette #11/12	SAGL11/12SS
Rigid Gracey Curette #13/14	SAGL13/14SS
Implant Curette (Universal)	SA4L/4RSST

Rigid Gracey Currettes*

SAGL5/6SS SAGL7/8SS SAGL11/12SS SAGL13/14SS

**Barnhart
Universal Curette**

SA1/2SS SA6/2SS

**Columbia
Universal Curette**


SA4L/4RSS SA13/14SS

**UC/Younger Good
Universal Curette**

SA7/8SS

**Loma Linda
Curette**

SAA10/11SS

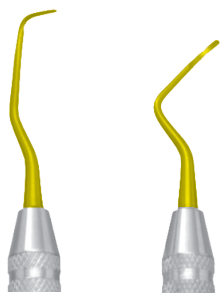Implant Currettes

SAGL11/12SST SA4L/4RSST

Sickle Scalers

SAMJ1/2NSS SA3/4SS SA6/7SS SA204SS

Explorers

XA11/12 XPW2F

Periodontal Probes

PAUNC12 P2NCP PACP12

CONE SOCKET SYSTEM**Cone Socket
Instruments****CARBON
STEEL** CSHA
3/8" CSW
Wrench**STAINLESS
STEEL** CSWS
Wrench**Universal Curettes**

R3 R4 4L 4R 7YG 8YG

R3S R4S 4LS 4RS 7YGS 8YGS

Sickle Scalers

S3 S4 S6 S7 MJ3 MJ4

S3S S4S S6S S7S MJ3S MJ4S

Rigid Gracey Currettes***CARBON** GL5 GL6 GL7 GL8 GL11 GL12 GL13 GL14 GL15 GL16**STAINLESS** GL5S GL6S GL7S GL8S GL11S GL12S GL13S GL14S GL15S GL16S

For a complete offering of
Hartzell instruments, call
1-800-950-2206
or visit www.ghartzellandson.com

*Also available in regular, extended and mini-extended

NOTE! If only one side of an instrument is shown, the
opposing side contains the same working end. Cone sockets
excluded, as they are supplied as individual working ends.

RESTORATIVE INSTRUMENTS







Restorative Instruments Premium Restorative Instruments

Hartzell's restorative instruments adhere strictly to Black's Operative Instruments Formulas. This comprehensive product line features premium non-stick composite instruments that will serve all of your restorative needs.



*Instruments are delivered pre-filled in the cassette and are not replaceable. Above photo is not representative of the assortment of instruments in the HZ-COMPKIT. Additional instruments to round out your armamentarium can be found on the opposing page.

Product	Item #
Composite Kit*	
Composite Kit*	HZ-COMPKIT
Contents*	
Cassette (10-Instrument)	IC10
Mirror	M5EA
Mirror Handle	CSH2A
Explorer	X5
Excavator	E38/39SS
Cord Packer	C113NS
Placement Instrument	CA6CH
Carver	CIPCT
Composite Instruments	P1/2T
	PUW/E4T,
	PWTNCIGFT3
Rubber Dam Forceps	IRDF

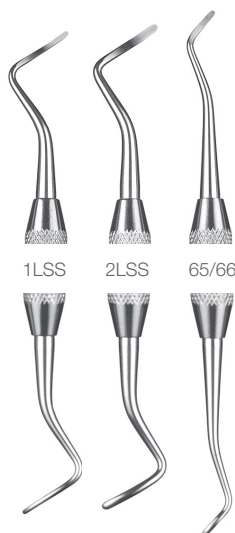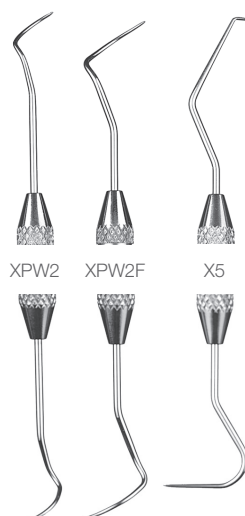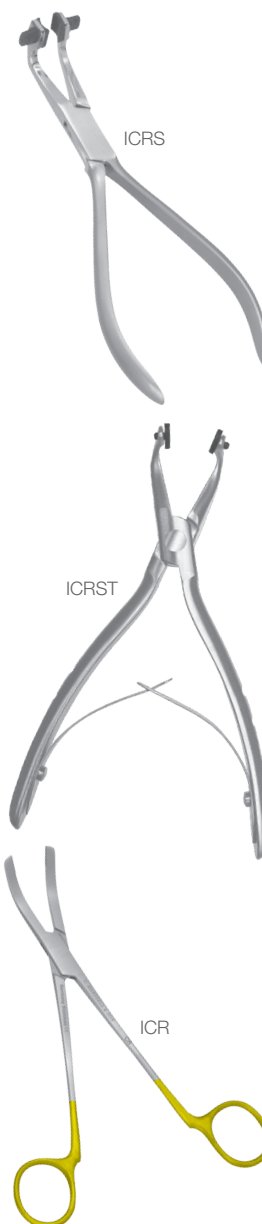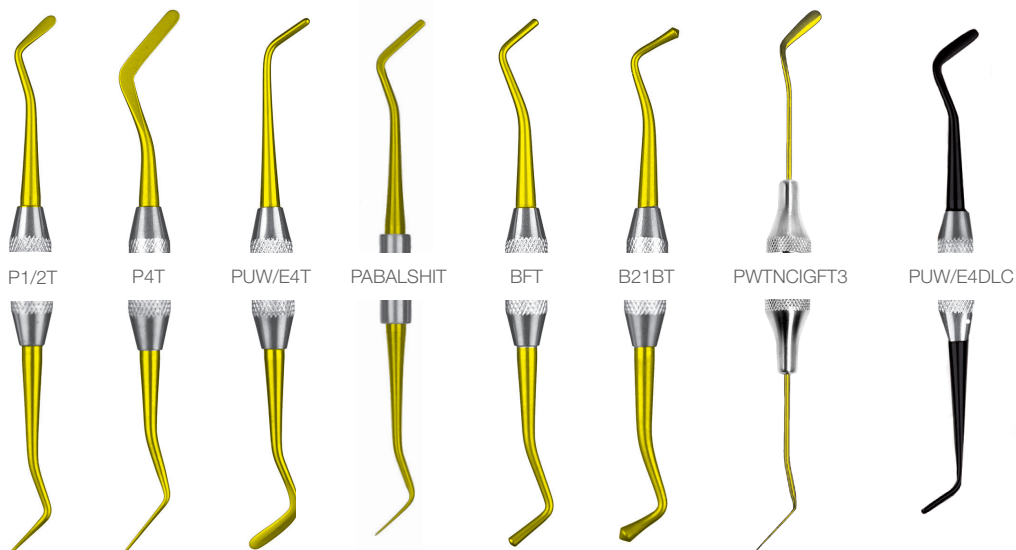**Spoon
Excavators****Excavators****Explorers****Crown Removers****Cord
Packers****Carver
Cord Packer****Carver****Hollenback
Carver****Interproximal
Carver****Discoid-
Cleoid****Ball
Burnishers****Composite
Instruments****Articulating Paper Forceps**

dm
DenMat®

SURGICAL INSTRUMENTS

Surgical Instruments

Top-Notch Surgical Instruments

Looking to get started with a set-up kit for implant dentistry? The recommended assortment below provides nearly everything you need to get started. Order part # HZ-IMPLANTKIT to order the entire kit, and save off of à la carte pricing! Additional instruments to round out your armamentarium can be found on the opposing page.



Periosteal — P9 (Molt)



Scalpel Handle — ISH5



Curette, Surgical Molt — SUM2/4 (Molt)



Bone Harvester/Scraper ISBS



Bone Condenser — CGAR3/5



Bone Condenser for Small Defects — C6/7




Probe, Periodontal (UNC 15mm) — PUNC15



Ruler — IRSS



Minnesota Cheek Retractor ICRM

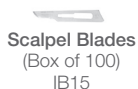

Scalpel Blades
(Box of 100)
IB15



Tissue Punch
6.5mm
IBP65L



Luxating Elevator
(2mm Curved)
IGLUX96



Bone
Mixing Bowl
IMBO



Irrigation
Solution Bowl
CUP6



Needle Holder
(Crile-Wood TC)
INCCW



Tissue Forceps
(Gerald Straight
Serrated)
ITPGS




Aspirator
(Frazier)
ASPF12



Suture Scissors
(Iris Curved)
ISIC



Ridge Mapping
(Bone) Instrument
IRMC

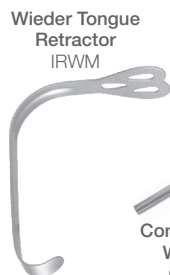

Wieder Tongue
Retractor
IRWM



Cone Socket
Wrench
CSWS



Lagrange Tissue
Scissors
ISL



Syringe,
Aspirating
SYRASPCP



Periotome


P1S

P2S


CSHA

ADDITIONAL INSTRUMENTS

Periotomes

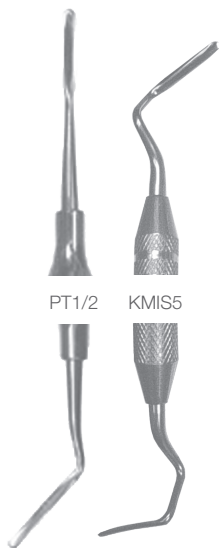

Periosteals

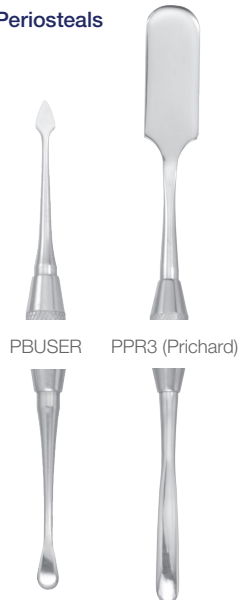

Surgical Curettes



Needle Holders (Tungsten Carbide)



Curved Scissors



Perio Knives



Elevators



Root Tip Picks




Luxating Elevators




Scalpel Handles



Extraction Forceps

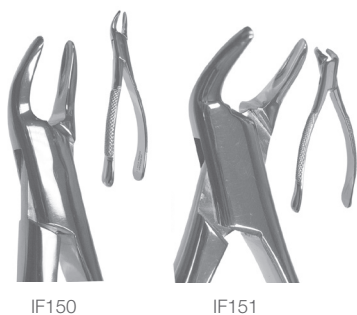

Deep-Gripping Extraction Forceps

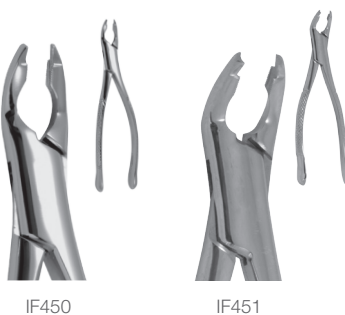


ENDODONTIC INSTRUMENTS

Endodontic Instruments

Precision Endodontic Instruments

Looking to get started with a set-up kit for endodontics? The recommended assortment below provides nearly everything you need to get started. Order part # HZ-ENDOKIT to order the entire kit, and save off of à la carte pricing! Additional instruments to round out your armamentarium can be found on the opposing page.

Hartzell
INSTRUMENTS

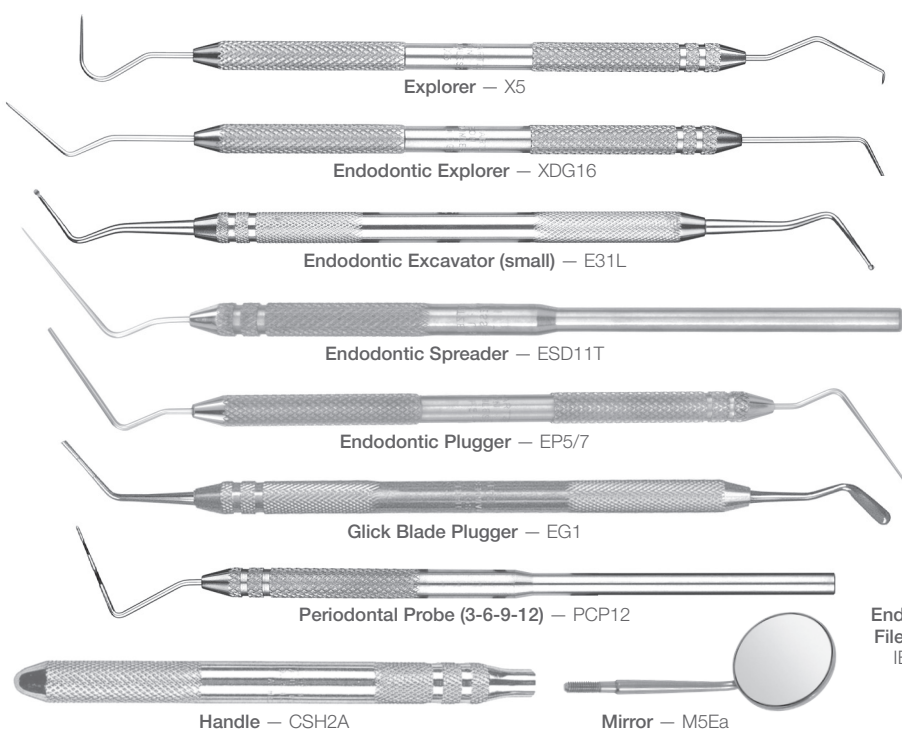


Rubber Dam
Forcep
IRDF



Endodontic Locking Pliers
(with grooves)
IPLE



Endodontic
File Holder
IEMFH

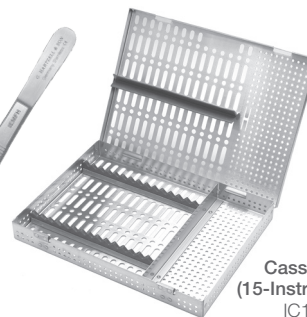

Cassette
(15-Instrument)
IC15



Rubber Dam
Clamp
IRDC1



Rubber Dam
Clamp
IRDC2



Rubber Dam
Clamp
IRDC9



Rubber Dam
Clamp
IRDC26

ADDITIONAL INSTRUMENTS


Endo Explorers



Endo Excavators



Endo Condensers

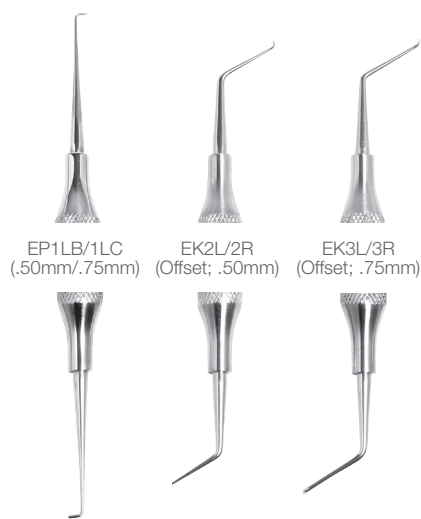

Periosteal



Dovgan MTA Carrier




Terauchi Gutta-Percha Removal Instruments



Dovgan MTA Pluggers



Perf Repair Instrument

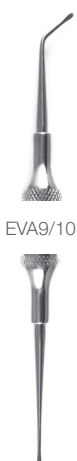

Huey Forceps (Diamond Coated) Available in 45° and 90° angles



Peets Forceps Available in 45° and 90° angles



Micro Mirrors



Simple Stem Handle



M.T.A. Pellet Forming Block

Patent #6,203,319
Designed by
Dr. Ed Lee,
San Francisco,
CA., U.S.A

CLEE
For transporting
and carving the
M.T.A. pellet.

IMP - Mixing Pads
for the Lee M.T.A.
Block are available
(50 pads).

For a complete offering of Hartzell instruments, call **1-800-950-2206**
or visit www.ghartzellandson.com



PeriOptix® lightweight frames and headlamps support better ergonomics, resulting in less strain on your neck and an overall more comfortable working experience.

PeriOptix is a leader in magnification and illumination innovations that include: a unique modular frame system, the first portable LED light source, a proprietary locking mechanism expressly designed to simplify the adjustment of flip-up loupes, and now the Firefly™ Cordless Headlamp System.

PeriOptix are available in through-the-lens and flip-up optics. We offer a full suite of high-quality magnification and illumination products, all at an unbeatable value.



NEW!

Firefly. The first and only cordless headlight designed to fit *your* loupes.

At only 28 grams, it is simply the most ergonomic headlight system available, offering superb brilliance, an optimized spot size, and extended run times.


FIREFLY™

3-position
adjustable temples

Adidas® Adivista frame
available in 8 colors

Rx-ready insert
available

2-Position
adjustable nose piece

3 types of optics
in 4 working
distances

 **PeriOptix®**

PeriOptix

Lightest Weight. Remarkable Value. Exceptional Quality.

Over 60% of Dental Practitioners suffer from back and neck pain, and over 90% suffer from eye strain. PeriOptix lightweight fitted loupes promote proper posture, reduce fatigue, enhance vision, and add to career longevity.



LED Headlamps

Lightest and Brightest Headlights on the Market

NEW!



Firefly Cordless Headlight System

The Cordless Headlight That Fits Any Pair of Loupes

- Unique cord-free design that eliminates the need for heavy power packs and annoying cords
- Allows you to easily upgrade your current loupes to incorporate a cordless headlight
- Easy to clean, reducing the risk of cross contamination
- Optimized spot size that eliminates the uncomfortable halo in your patients' eyes while providing uniform and natural-color output

Product	Item #
Firefly Cordless Headlight System	
Firefly Cordless Headlight System	LED5
<i>Contents: 2 rechargeable light pods, magnetic loupe mount and hardware, dual charging station for light pods, cleaning cloth, hook-resistant carrying case, instructions</i>	
MT-S Flip-Up Light Mounts	
Panoramic	MT-P3
Legacy	MT-P2
Orasoptic® Square (old)	MT-O
Orasoptic Round (new)	MT-O-N
Sheer Vision	MT-P-O
Surgitel®	MT-P3
Zeiss®	MT-Z
TTL Light Mounts	
Adivista	MT-A
Small Universal	MT-U-S
Large Universal	MT-U-L
Orasoptic Titanium Frame	MT-Ti
Small TTL Solid	MT-T-S3

MicroLine™ Mini LED Headlight

The Lightest LED Headlight

- Weighs less than nine grams
- Brightest mini LED headlight (40,000 Lux*)
- 8-hour+ rechargeable, replaceable Li-Ion battery
- Attaches to all PeriOptix and most competitor loupes or Freedom™ headbands
- 90% charge in only 90 minutes
- User-replaceable cable
- Composite filter available



MicroLine Mini on MicroLine TTL

Optional composite filter



Only PeriOptix lights have an inexpensive, replaceable cable

Solaris® Portable LED Headlight

Brightest Lightweight LED Headlight

- Brightest LED headlight (80,000 Lux*)
- 5-hour+ rechargeable, replaceable Li-Ion battery
- Attaches to all PeriOptix and most competitor loupes or Freedom headbands
- 90% charge in only 90 minutes
- User-replaceable cable
- Composite filter available

Solaris on Hogies® System



PeriOptix power pack features a snap-in battery

*Lamp Illumination 80,000 lux @ 35 cm (7,500 foot-candles @ 14"), 40,000 lux @ 35 cm (3,750 foot-candles @ 14")

To order, call **1-800-4DENMAT** (1-800-433-6628) or visit www.denmat.com

dm
DenMat®


Certified PeriOptix® Loupe Specialist

A Certified PeriOptix Loupe Specialist can help you get fitted with the perfect flip-up or through-the-lens loupes



Pupillary Distance


Pupillary Distance (PD) is the measurement, in millimeters, between the centers of the pupils in each eye. This measurement can be done quickly and painlessly by your local PeriOptix representative or by any optometrist. Note: If measured by an optometrist make sure the measurement is taken at infinity (∞).

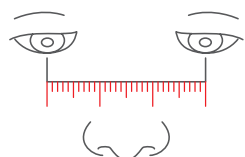

Right PD	Total PD	Left PD

Working Distance


Working Distance is the measurement at which the user is comfortable working. This distance should be calculated while the doctor is in an upright, and comfortable position that keeps the clinician's ergonomics in check.

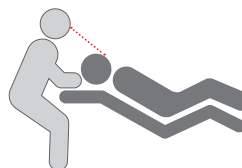


S	12 – 16 in (30 – 40 cm)
R	14 – 18 in (35 – 45 cm)
L	16 – 22 in (40 – 55 cm)
XL	18 – 22 in (45 – 55 cm)

Field of View for PeriOptix A reference for magnification and field of view


1.0x

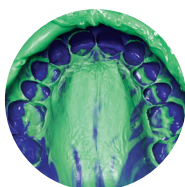

2.5x

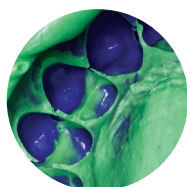


3.5x

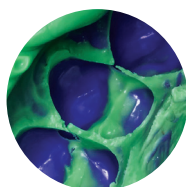


As seen on a SplashMax impression.

Small Field of View

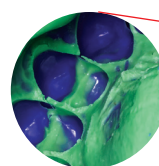

Large Field of View

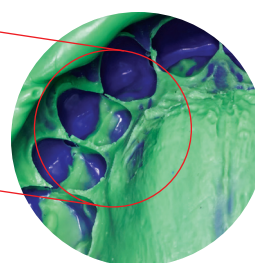

Higher magnification typically results in a smaller field of view.

Through-the-Lens (TTLs) Loupes

Premium, Lightest-Weight TTL Loupes on the Market



2.3x MicroLine™ TTL Loupes

The Lightest Option on the Market

- Extremely lightweight at only 36 grams
- Lightest and smallest loupes available
- Custom-fit to your specifications



PeriOptix 2.3x MicroLine optics are less than 1/3 the size of conventional optics
33% SMALLER!



2.7x TTL Loupes

A Perfect Balance of Clarity, Magnification, and Field of View

- Seven standard interpupillary distances
- Optimum clarity and depth of field



PeriOptix 2.7x TTL optics are much lighter than conventional optics



3.1x TTL Loupes

High-Magnification TTL

- Custom-fit high-magnification optics
- Highest declination angle



PeriOptix 3.1x MicroLine optics offer the greatest level of magnification

Flip-Up Loupe Optics

Lightweight and Versatile Loupes

Panoramic Flip-Up Loupes

Our Highest-Rated Flip-Up Loupes

- Three magnification options: 2.5x, 3.0x, and 3.5x
- New streamlined design
- Panoramic wide field lenses
- Waterproof optics
- Rx lenses available

Approximate Magnification	Size	Approximate Working Distance	Field of View
2.5x	S	11" – 14" (28 – 36 cm)	3.5" (8.8 cm)
	R	13" – 17" (33 – 43 cm)	4.1" (10.5 cm)
	L	15" – 20" (38 – 51 cm)	5.1" (12.9 cm)
	XL	17" – 22" (43 – 56 cm)	5.7" (14.5 cm)
3.0x	S	11" – 14" (28 – 36 cm)	2.7" (6.8 cm)
	R	13" – 17" (33 – 43 cm)	3.5" (8.8 cm)
	L	17" – 21" (43 – 53 cm)	4.1" (10.5 cm)
	XL	19" – 24" (48 – 61 cm)	2.4" (6.1 cm)
3.5x	R	14" – 17" (35 – 43 cm)	3.5" (8.8 cm)
	L	17" – 21" (43 – 53 cm)	2.9" (7.3 cm)
	XL	19" – 23" (48 – 58 cm)	3.1" (7.8 cm)



Legacy™ Sport Flip-Up Loupes

Greatest Value on the Market

- Three magnification options: 2.5x, 3.0x, and 3.5x
- High quality value line loupes
- Multiple frame choices



MicroLine™ Flip-Up Loupes

Extremely Lightweight with Exceptional Clarity

- Two magnification options: 2.3x and 2.7x
- Fully adjustable
- Lightest weight optics
- Rx lens inserts
- Large/Wide field of view

Approximate Magnification	Size	Approximate Working Distance	Field of View
2.3x	S	11" – 15" (28 – 38 cm)	3.3" (8.3 cm)
	R	14" – 18" (35 – 45 cm)	4.0" (10.3 cm)
	L	16" – 22" (40 – 55 cm)	4.8" (12.2 cm)
	XL	19" – 25" (47 – 63 cm)	5.6" (14.1 cm)
2.7x	S	11" – 15" (28 – 38 cm)	2.2" (5.7 cm)
	R	14" – 18" (35 – 45 cm)	2.8" (7.0 cm)
	L	16" – 22" (40 – 55 cm)	3.3" (8.3 cm)
	XL	19" – 25" (47 – 63 cm)	3.8" (9.6 cm)



To order, call **1-800-4DENMAT** (1-800-433-6628) or visit www.denmat.com

Frames

Light, Comfortable and Stylish Frames

PeriVista Frame by Adidas®

Our Lightest, Most Comfortable Frames

- Lightweight Adidas Adivista sports frame
- Three-way adjustable temples
- Two-position nosepiece
- Four working distances
- Available in eight distinct frame colors



Matte Black



Red



Grey



Copper



Gold (Titan)



White



Blue (Petrol)



Silver



Laser Inserts
Now Available



Legacy Sport Frames

Our Greatest Value

- Multiple frame colors
- Rx option (Mirage only)



Mirage
with available Rx insert
and Legacy flip-up loupes



Red Mirage



Black Mirage




Red Rimz



Black Rimz




Pink Rimz

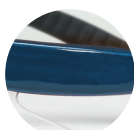


Blue Rimz



White Rimz

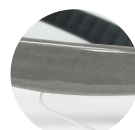

Silver Rimz

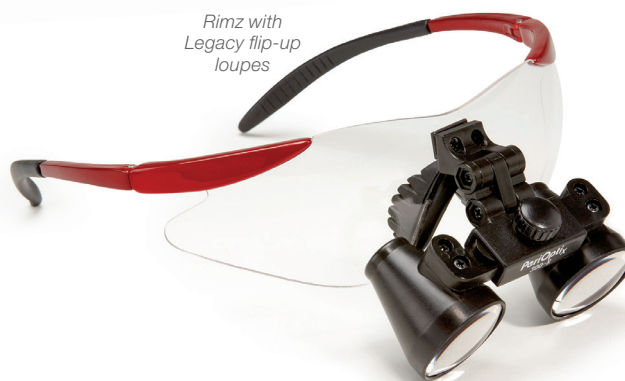

Rimz with
Legacy flip-up
loupes

Flip-Up Loupes/Frame Systems

Lightweight and Versatile

Air-X Flip-Up Loupes

Lightweight Loupes Ideal for Laser Dentistry

- Can be worn over eyewear
- Snap-on safety shield
- Adjustable nose pad
- Optional laser shield inserts available



Includes snap-on shield for complete splash and U.V. protection. Shown with Legacy Flip-up Loupes.




Air-X is designed to be comfortably worn over existing eyewear

Laser Source	Wavelength (Nm)	Green (F-U-A-Sg)	Pink (F-U-A-Sp)
Diode	808/810, 940, 980, 1064	•	•
Nd:YAG	1064	•	
CO ₂	10600	•	



Laser inserts are ideal for use with the NV Diode Laser from DenMat.

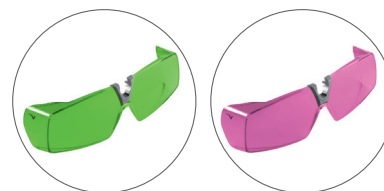

Optional laser shield inserts

Hogies® Advanced Optical Platform

Extremely Versatile, Magnetic Frame System

- Can be worn over eyewear
- Magnetic modular sports frame system
- Prescription lens option
- Replaceable shields and nose pads




Hogies with magnetic loupe attachment and Legacy Flip-up Loupes




Hogies are designed to be comfortably worn over existing eyewear



Red

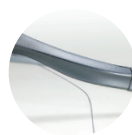

Smoke Grey

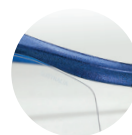

Metallic Blue




Black

To order, call **1-800-4DENMAT** (1-800-433-6628) or visit www.denmat.com

dm
DenMat®



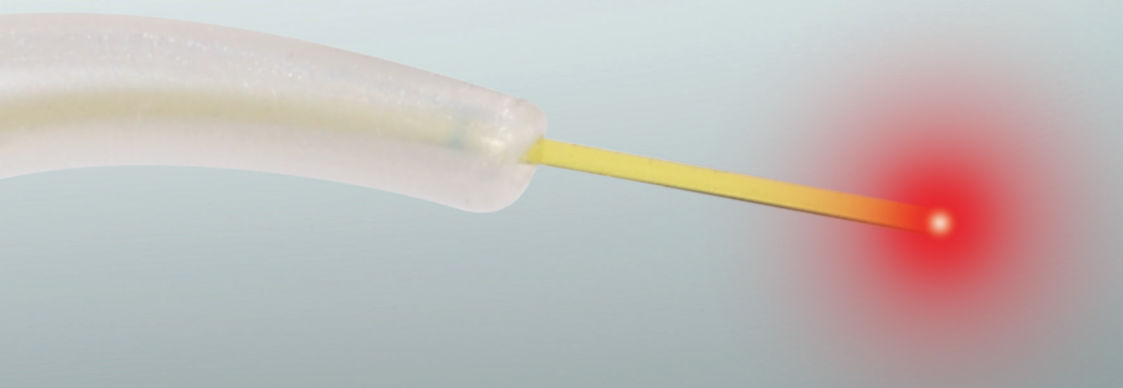

From our revolutionary NV[®] Microlaser to the best-in-class Sapphire[®] Curing Light, DenMat[®] offers a wide array of award-winning equipment to help build your practice.



NV Microlaser

Cordless Soft-Tissue Diode Laser

The NV Microlaser delivers the power and capabilities of larger desktop soft-tissue lasers in a pen-sized package. Weighing only 1.9 ounces and completely self-contained, NV diode laser's wireless hand-held design gives you unprecedented freedom in the operatory. Experience the benefits of soft-tissue laser technology without any additional clutter or cords. Experience NV Microlaser.





Product	Item #
NV Microlaser	
NV Microlaser	ZLR1001
NV Microlaser Accessories & Replacement Parts	
NV Charger Base	ZLR1006
NV Wireless Foot Pedal	ZLR1008
NV Microlaser Lithium Ion Battery	ZLR1010
NV Disposable Regular Tips 5mm (25-pack)	ZLR1012A
NV Disposable Perio Tips 7mm (25-pack)	ZLR1013A
NV Wireless Power Cord North America (wall plug only)	ZLR4024
NV Microlaser Power Supply (does not include wall plug)	ZLR1037
Universal Diode Laser Safety Glasses	LR1008

Looking to perform more predictable pain-free dentistry?

DenMat's soft-tissue diode lasers promote safer, healthier tissue management for your patients while increasing efficiency in your practice.

Benefits of Laser Dentistry Include:

- Minimized need for topical anesthetic
- Improved access and field of vision
- Minimal tissue damage
- Minimal post-op sensitivity
- Reduced healing time and treatment time in many cases
- Keep procedures in house and enjoy better patient compliance
- Incremental revenue through Laser Assisted Periodontal Therapy (LAPT)
- Plus many more!


Frenectomy



Gingivectomy



Perio (LAPT)

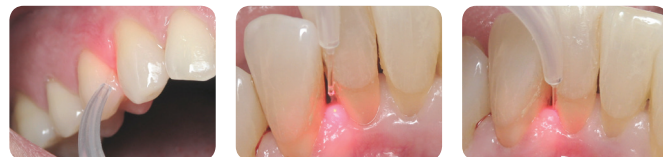

Laser Troughing





SOL® Desktop Laser

Portable Soft-Tissue Diode Laser

SOL delivers true portability, power and precision, along with enhanced ergonomics, simple interface, and a high-contrast aiming beam designed to make laser dentistry easier than ever. It also happens to be incredibly affordable. Great curves. Outstanding performance. Unbeatable value.

- Easy-to-see blue aiming beam
- Easy-to-use disposable tips
- Max power: 3 Watts in continuous mode, 5 Watts in pulsed mode

Product	Item #
SOL Portable Diode Laser	
SOL Laser	LR1001W
SOL Laser Power Cord	LR1005
SOL Replacement Fiber Optic Cable	LR10006
SOL Wireless Foot Pedal Upgrade Kit	LR1045
Diode Laser Tip 4mm, 400µm (25-pack)	033980028-0
Diode Perio Tip 9mm, 300µm (25-pack)	LR1010
Universal Diode Laser Safety Glasses	LR1008


Visibility




True Portability



Simplicity

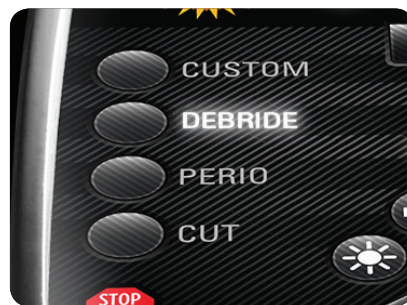


Wireless Foot Pedal

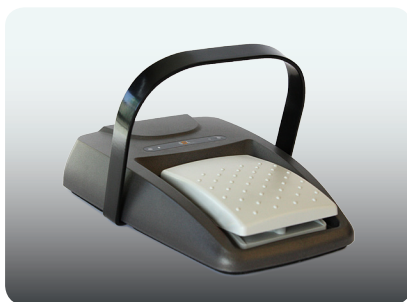

Ergonomics



Single-Use Tips

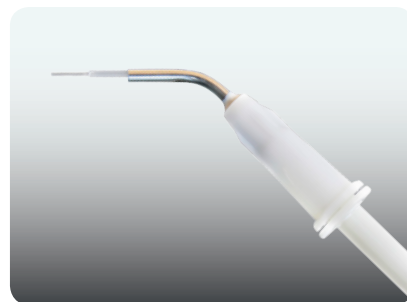



Sapphire® PAC Curing Light

The Sapphire light is a versatile and powerful plasma-arc curing and whitening system, all in one package. Light-cured composites cure in 3 – 5 seconds and dual arch bleaching is completed in one office visit. Fast, convenient, comfortable, economical... unparalleled.

- No harmful UV rays and no need for skin protection
- Whitens individual teeth as well as full arches with the Sapphire Whitening Crystal
- Super-fast curing in 3 – 5 seconds, saving an estimated \$40,000/year in chairtime²
- Curing power that goes the distance. Effective curing from up to 10mm with greater power than any other light tested from any distance¹
- Fully cures resins in less time than our competitors¹
- Reduces finishing time for Lumineers up to 50%

Product	Item #
Sapphire Plus	
Sapphire Plus Light (US)	033968000
Sapphire Plus Light (Canada)	033968059
Sapphire Components	
Sapphire 2mm Ceri-Taper Tip	033958506
Sapphire 4mm Ceri-Taper Tip	033958501
Sapphire 9mm Tip	033958519
Sapphire 12mm Tip	033958503
Whitening Crystal	030395250
Sapphire Wall Mount Bracket	043958417
Rolling Stand w/Arm	030597275
Rolling Stand w/No Arm	043958410
SBP Arm Only with 360-Degree Articulation	030595375
Sapphire Power Cord	043958045
LT Scope Hand Piece Holder Assembly	044391059
Sapphire Plus Light Guide with Pistol	033968553
Sapphire Plus Replacement Handle – Right Side	043960937
Sapphire Plus Replacement Handle – Left Side	043960936
Sapphire Tip Receptacle	033958080
LT Orange Shield w/Grommet	043957515
Sapphire Plus Extended Warranty	900039900
Sapphire Lamp Module	043951176
Keys	033968859



Flashlite Magna® 4.0 LED Curing Light with Advanced Beam Collimation

Flashlite Magna 4.0 features an LED engine that delivers optimal beam collimation, reducing divergence and maximizing curing power in deep preps.

- 1300 mW/cm²
- 10mm beam width
- Short head angle allows easy access to posterior teeth
- Ideal for large restorations
- Cordless
- 2-year warranty

Product	Item #
Flashlite Curing Lights	
Flashlite Magna 4.0	CR1079
Flashlite 2.0	CR1080
<i>Contents: curing light, AC/DC wall plug adapter, light shield (4.0 kit only), replacement lens caps (3pc), tacking tips (5pc), charger base with built-in radiometer</i>	
Flashlite Accessories & Replacement Parts	
Charger Base for Flashlite Magna 4.0	CR1083
Charger Base for Flashlite 2.0	CR1063
Tacking Tips (5-pack)	CR1075
Replacement Lens Caps for Flashlite Magna 4.0 (25-pack)	CR1082
Replacement Lens Caps for Flashlite 2.0 (25-pack)	CR1024
Light Shield	CR1077
Barrier Sleeves (500-pack)	CR1078



Flashlite® 2.0 LED Curing Light

The cordless, ultra compact design, exceptional quality and value price make the Flashlite 2.0 an ideal LED curing light for every operator.

- 1100 mW/cm²
- 7.2mm beam width
- Short head angle allows easy access to posterior teeth
- Ideal for curing small restorations and sealants
- A perfect hygiene light
- 2-year warranty



The NEW DenMat Lab will partner with you to create one-of-a-kind smiles using the latest minimally-invasive technology and techniques.

Send us your challenging cases and let us guide you with our new A.R.T. process, which helps make the complex simple. Signature Service is our new premium level of esthetics designed to meet your high standards.

Lumineers® and Snap-On Smile® are two of our most popular brands. Our experienced technical consultants are here to help you with design, product information and specifications, as well as any other questions you may have.

Let us partner with you for all your esthetic needs.

LUMINEERS®

The #1 patient-requested veneer

**Veneers**

DenMat Lab, the innovator of thin veneers with over 30 years of experience, is here to help you with your challenging esthetic cases. Lumineers technology and processes allows you to provide beautiful smiles with the greatest preservation of tooth structure possible. We now offer our new A.R.T. and Signature Service processes, as shown on the next page. Send us a case and allow us to show you why Lumineers is the #1 patient requested veneer.

Materials

- Cerinate II lucite-reinforced porcelain
- Lithium disilicate

Preparation Techniques

- No prep
- Minimal prep
- Traditional prep

Fabrication Techniques

- Pressed
- Layered

LumiZir®

Solid Beauty

**Permanent Restorations**

LumiZir is our brand of zirconia, available in full-contour restorations made from one monolithic piece of zirconia. LumiZir Natural offers high translucency and 1160 MPa strength. We also offer LumiZir stacked with porcelain and lithium disilicate for a more customized esthetic case.

We offer several materials for permanent restorations that will allow you to customize your cases to optimize beauty and strength.

Now available in:**Crowns**

- LumiZir/LumiZir Natural™
- Lithium disilicate
- Porcelain fused to metal (PFM)
- Captek
- Full Cast

Bridges

- LumiZir
- Lithium disilicate
- Porcelain fused to metal (PFM)
- Maryland

Inlays/Onlays:

- Lithium disilicate

snapon
smile®

Changing Smiles...Changing Lives.

**Provisionals**

Many dentists regard temporization as one of the most difficult aspects of full-mouth restoration. Snap-On Smile is not only an excellent provisional, it also provides patients with exceptional esthetics. It can last for years, yet is affordable enough to be temporary. For patients who are hesitant to commit to more involved restoration treatments, Snap-On Smile gives them the experience of living with better function and superior esthetics.

- Full or partial arch
- Cosmetic enhancement
- Implant provisional
- Removable partial denture
- Non-invasive

97% Excellent or Good by Clinicians Report†

81% of 31 CR Evaluators stated they would incorporate Snap-On Smile into their practice.

Snap-On Smile was primarily used by CR Evaluators for:

- Provisional restoration
- Short-term esthetic improvement
- Long-term esthetic improvement
- Diagnostic set-up

*Data on file.

†An independent, non-profit dental education and product testing foundation: Clinicians Report,® February, 2012. For the full report, contact DenMat.



A.R.T. (Additive-Reductive Template) is DenMat's new, innovative process designed to make even complex veneer cases simple. This revolutionary process provides a "blueprint" from the Lab to the dentist to ensure optimal, easy to follow, conservative prepping for each individual case.

A.R.T. is available in 3 options to suit all: no prep, minimal prep, and idealized prep. Idealized prep is for advanced cases or combo cases that include crowns, bridges, implants, etc.



Signature Service is a premium level of esthetics and service specifically designed to provide superior results, even for the most demanding cases.

- Case supervision by master ceramist Brad Jones, FAACD
- Advanced color-matching
- Enhanced artistry and nuances
- State-of-the-art processes for working with porcelain
- Smile-design technical specialists with advanced training

RSVP Service

- Reserve your case in our production schedule
- Request a specific technician to design all of your cases with a reservation
- No-charge digital photographs during wax or completed phase
- Able to apply various occlusal philosophies to case design (Dawson, LVI, OBI, and Pankey)
- Lumineers with Signature Service Pricing: \$279 per unit vs. \$239 per unit for standard Lumineers – boutique lab results for a better value!



Example of an Advanced A.R.T. Case



Pre-op Model



Facial Reduction



Interproximal Reduction



Incisal Reduction



Gingival Reduction



Additive/Reductive Wax-up



Matrix of Wax-up for the "Trial Smile"

Product	Item #
Kits & Replacements	
Lumineers Placement Kit	Lum1005
Snap-On Smile Starter Kit	056406000
Snap-On Smile Carrying Case	056404101B
Snap-On Smile Adjustment Bur Kit	033898100
Snap-On Smile Cleaning Solution	056404004
Snap-On Smile Anti-Bacterial Gel	056404006



Lumineers Placement Kit



Snap-On Smile Adjustment Kit



Snap-On Smile Starter Kit



Snap-On Smile Carrying Case



Snap-On Smile Cleaning Solution



Snap-On Smile Anti-Bacterial Gel

Product	Item #
Practice Enhancement	
Lumineers Patient Model	050405175
Photography Kit	040750500

Professional Dental Photography Kit



To Get Started

Call 1-800-433-6628 ext. 5520 and ask customer service for prescription forms, shipping labels and boxes or ask our technical lab experts for a free consultation.

To order, call **1-800-4DENMAT (1-800-433-6628)** or visit www.denmat.com





LumiGrip®


Veneer Placement Tool

The gentle way to place Lumineers, or any veneer, with ease and confidence! The sleek LumiGrip is an ergonomically designed, fully autoclavable, anodized aluminum handpiece that maximizes suction forces from standard saliva ejector and HVE lines in order to securely and delicately assist in the stress-free and precise placement of “contact lens-thin” Lumineers. Held like a pen, the powerful suction force directed through the soft, baffled, silicone cup securely grips the veneer, enabling precise placement. An effortless and slight retraction of the index finger against the rotary valve on the body of the handpiece is all that is needed to turn off the suction and release the veneers.

LumiGrip Kit



LumiGrip Handpiece

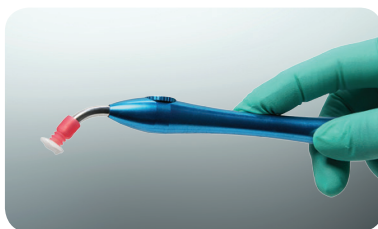


Product	Item #
LumiGrip Kit	
LumiGrip Kit	031752150
<i>Contents: LumiGrip handpiece, LumiGrip suction tips (25pc), LumiGrip cleaning brush, LumiGrip grey hose</i>	
LumiGrip Accessories	
LumiGrip Suction Tips (10pc)	030752110
LumiGrip Grey Hose	040752102



Ceri-Saw & Ceri-Sander


Interproximal Finishing Tools

These tools are particularly useful for interproximal finishing of multiple veneers and separating fused composite restorations. The ideal handle allows for easy manipulation of the thin blades that give you total control to prevent accidental slips and injury. Designed for fast, easy-to-use, interproximal finishing, the Ceri-Saw and Ceri-Sander are safe and more comfortable for the patient than traditional disc and strip systems.

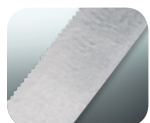

Ceri-Saw

The Ceri-Saw is used for interproximal finishing of multiple laminates and separating fused composite restorations. Its exceptionally thin, sharp saw blades can easily cut through composite without concern for cutting porcelain.

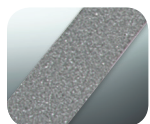

Ceri-Sander

The Ceri-Sander is a single-sided, diamond coated strip in fine and ultra-fine textures. The Ceri-Sander is for final finishing and smoothing of the interproximal.



Product	Item #
Ceri-Saw Anterior Kits & Blades	
Ceri-Saw Anterior w/11 Blades <i>Contents: anterior handle w/11 blades</i>	031336500
Ceri-Saw Anterior w/1 Blade <i>Contents: anterior handle w/1 blade</i>	031336510
Ceri-Saw Anterior Blades (10pc)	031336520
Ceri-Saw Posterior Kits & Blades	
Ceri-Saw Posterior w/11 Blades <i>Contents: posterior handle w/11 blades</i>	031336600
Ceri-Saw Posterior w/1 Blade <i>Contents: posterior handle w/1 blade</i>	031336610
Ceri-Saw Posterior Blades (10pc)	13371
Ceri-Sander Anterior Kits & Blades	
Ceri-Sander Anterior Kit <i>Contents: anterior handle, fine strips (6pc), ultra-fine strips (6pc)</i>	031337500
Ceri-Sander Fine Strips Large (6pc)	031337510
Ceri-Sander Ultra-Fine Strips Large (6pc)	031337520
Ceri-Sander Posterior Kits & Blades	
Ceri-Sander Posterior Kit <i>Contents: posterior handle, fine strips (6pc), ultra-fine strips (6pc)</i>	031347500
Ceri-Sander Fine Posterior Strips Small (6pc)	031347510
Ceri-Sander Ultra-Fine Posterior Strips Small (6pc)	031347520



Lumineers® Finishing Kit

Veneer Placement Finishing System

This finishing system is specially designed for Lumineers placement. However, it's a great set of tools for any veneer placement. This all-inclusive kit makes the finishing process quicker and much smoother than do other tools on the market today. Each tool is specifically designed to finish the laminate surface and interproximals without causing the patient pain or damaging the laminates.



Matrix Strips

Interproximal Strips

Matrix Strips enclose the restoration, creating a temporary interproximal surface to prevent the spread of filling material.

Product	Item #
Lumineers Finishing Kit	
Lumineers Finishing Kit	033897100
<i>Contents: Lumineers Diamond and Bur Block, Schure 349 Instrument, porcelain laminate polishing paste syringe 5g, Ceri-Saw and Ceri-Sander handles (1ea), Ceri-Saw Anterior blades (10pc), Ceri-Sander fine strips (6pc), Ceri-Sander ultra-fine strips (6pc), syringe tips 18G (10pc)</i>	

Product	Item #
Matrix Strips	
Matrix Strips (120pc)	031000300



Thriving In This Economic Environment



Our Economy



Your

At DenMat®, we firmly believe that we are only successful if we make you, our customers, more successful practitioners and business owners. To that end, we offer a full suite of practice-building and continuing education services designed to elevate your standard of care and grow your practice.

The Lumineers Smile Discovery Program® (LSDP) is designed to ensure that your practice can take maximum advantage of the extensive consumer brand awareness of Lumineers. From targeted national advertising campaigns to customized marketing support and tools to assist in case acceptance, LSDP is your machine for practice growth.

DenMat Education provides a wide variety of innovative and comprehensive education programs for dental practitioners. These educational programs include national conferences, regional hands-on training courses, live patient placement courses, webinars and more. DenMat Education strives to provide educational environments that support truly practice-changing learning that not only enhances clinical practice, but also has a powerful impact on a practice's bottom line.



LSDP — The Premier Practice-Building Program in Dentistry

The Lumineers Smile Discovery Program (LSDP) is designed to ensure that your practice can take maximum advantage of our consumer brand awareness. From targeted national TV advertising campaigns to customized print ads and emails, LSDP is your machine for practice growth. You'll get what you need to drive your practice to new levels, including:

- Tools for case acceptance
- Patient demand generation
- Comprehensive practice marketing support
- FREE monthly product



LumiSmile® — The Key to Case Acceptance

Give your patients the unique experience of seeing their new, perfect smile — instantly — with the LumiSmile Program. Your patients are up to 500% more likely to say yes to Lumineers or other cosmetic procedures after seeing their LumiSmile Digital Smile Makeover.



National Advertising

Targeted national TV advertising and on-line pay-per-click marketing are designed with one objective in mind: drive new patients to your office. The campaign drives literally millions of impressions a year through consumer media. In fact, Lumineers.com is often the #1 ranked site for "veneers," the top dental search term on Google. Our website and doctor locator capture interest and make sure new patients discover your practice.

"It's like we've got a national advertising agency working for our practice!"

— Dr. Peter Harnois, Hinsdale, IL



LSDP Member's Website

The LSDP Portal gives you access to easy-to-use ad templates and emails that you can customize for your practice, phone scripts for your team to use to communicate your services and practice-growing marketing tips and tools.

Trained Lumineers Smile Discovery Program consultants will help set up and train your office to make sure you're getting the most out of this valuable aspect of the program.



LumiDays

LumiDays are sponsored Open House events held in your practice to show potential patients what you and Lumineers can do for their smiles. We support your LumiDays with promotional flyers and actual LumiSmile makeovers during the event. It's a proven way to drive new cases!



You'll receive FREE monthly product as a member of LSDP _____

LumiSmile
White

LUMI BRITE

pro-denRx

ATRIDOX ViziLite

Splash!

Splash!
MAX

precision

Bite
Registration

PERFECTray

The GRIPPER

perfec
temp10
Temporary Crown & Bridge Material

PERFECTEMP II

nuance
COMPOSITE DISPENSER

cp
Core Paste

Infinity

Tenure

Ub
Ultra-Bond

Virtuosos

Geristore

dm
DenMat®

To order, call 1-800-4DENMAT (1-800-433-6628) or visit www.denmat.com

Continuing Education

Whether you want to increase your esthetic production, hygiene production, or both, we have a course that's right for you. Learn what's new *and* how to incorporate it into your practice — starting Monday morning!

Soft Tissue Management®

A healthy smile is your business...
helping you succeed is ours.



Would you like to:

- 1 Get more patients to say "YES"?
- 2 Increase hygiene production while seeing fewer patients?
- 3 Empower your hygiene team to be the driver of elective procedures?
- 4 Adapt the latest technologies and techniques for non-surgical perio therapy?
- 5 Energize, motivate and inspire your entire office staff?

Then add one of DenMat's Soft Tissue Management (STM®) Continuing Education courses to your calendar. This 1-day workshop features tried and proven methodologies that will revolutionize your practice and turn your current hygiene program into a profit center. The latest technologies and techniques will be presented that will teach your staff how to perform non-surgical perio therapy and bill appropriately for the services you are already providing to your patients. You will learn first-hand from one of a select group of speakers that

have successfully implemented this program in their own practices, and leave the course with a new perspective and motivation to make Soft Tissue Management a more central part of your service offering.

Your patients will thank you. Your team will thank you. And your hygiene production will transform. Get inspired at a DenMat STM course near you.




Experience one-on-one demos of STM products and techniques!



Enhanced visualization with
PeriOptix® Loupes & Lights



Heated sub-gingival irrigation
with the PSP® Scaler

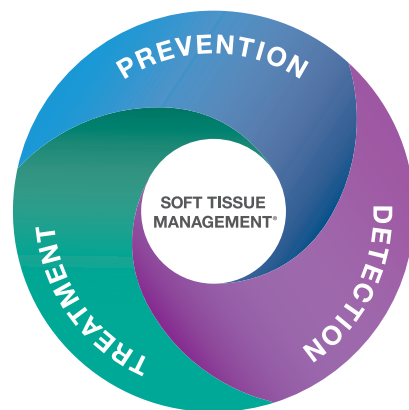

To order, call **1-800-4DENMAT** (1-800-433-6628) or visit www.denmat.com

Mastering Minimally-Invasive Smile Design



Join us for two days of powerful continuing education and learn:

- Best practices from your peers who have built highly successful practices with minimally-invasive esthetics in a challenging economy.
- Advanced esthetic dentistry, including case selection, smile design, preparation using the new Additive-Reductive Template (A.R.T.), and final finishing and seating of minimally-invasive veneers.
- Practice marketing, patient communication, case acceptance and practice management skills to build a successful esthetics-oriented dental practice.

13
CE CREDITS



Before




After



Before



After

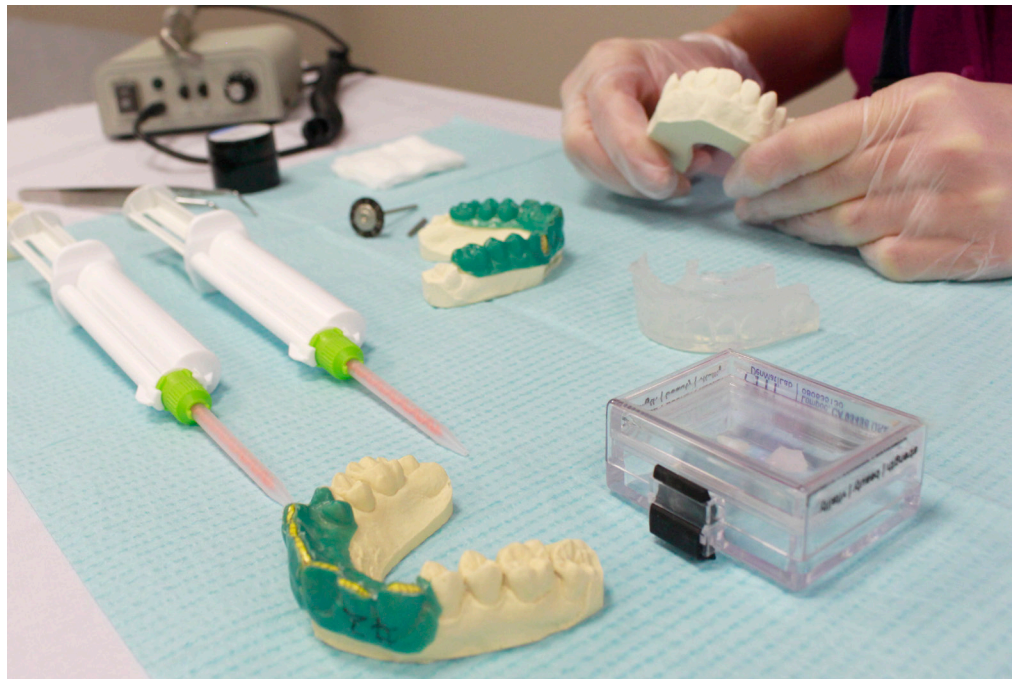

Learning Objectives

- Learn both clinical and business implementation all in one course – all you need to make a difference in your practice right away!
- Experience engaging and unique hands-on training with different veneer case types and small equipment and material demonstrations.
- Discover the new additive-reductive template (A.R.T.) and process to minimize chair time and enable the most conservative preparation with the best esthetic outcome.
- Make complex cases simple — all with minimal tooth reduction. Learn to tackle challenging alignment and asymmetrical arch forms, and many other complex cases.
- Participate in class discussions about case planning; bring your cases for personalized evaluation!



DenMat Education is designated as an Approved PACE Program Provider by the Academy of General Dentistry. The formal continuing education programs of this program provider are accepted by AGD for Fellowship, Mastership and membership maintenance credit. Approval does not imply acceptance by a state or provincial board of dentistry or AGD endorsement. The current term of approval extends from 8/1/2014 to 7/31/2017. AGD Provider #218048

For more information and a full list of DenMat Education courses, please visit www.denmat.com/courses.

LUMINEERS®

snayon
smile

dm
DenMat® EDUCATION

To order, call **1-800-4DENMAT** (1-800-433-6628) or visit www.denmat.com

On-line shopping with DenMat

Enjoy the Benefits of On-Line Ordering

DenMat customers who order on-line receive these benefits:

- Save Money – Registered users receive available promotion codes and free shipping on certain orders.
- Save Time – Ability to view and select “Your last order.”
- Availability – Website is open 24/7, 365 days a year.
- Convenience – “Auto-filled” forms for logged-in registered users, “Shopping List” to save selected products for later purchase, available saved multiple selection “Ship to Addresses.”
- Freedom – Choose your payment type and shipping.
- Security – Encrypted secured data protection keeps your data safe and secure.

DenMat's On-Line E-Commerce Guide

Follow the easy registration steps on the following page to discover the many benefits of ordering on-line with DenMat's e-commerce portal.

- How to “Sign up”
- How to “Log in”
- How to “Place an Order”
- Benefits of “On-line Ordering”
- Product Search Feature
- Customer Service
- E-Commerce Demo Video

SL3 Diode Laser

Welcome, DENMAT DENTIST | Log Out

Forgot Password? | New Customer?

1-800-4DENMAT (1-800-433-6628) | 1-805-346-3700

Downloadable Forms

Products DenMat Lab Marketing DenMat Education Portal Fast Order

Products : Soft Tissue Lasers : SL3 : SL3 Diode Laser

Print this page Email a Friend

SL3® Diode Laser

Overview Order Now Videos Product Specifications

Part#	Product Description	Qty	My Price	LSDP Price	
ZLR2002	SL3 Desktop Laser	1	\$3,999.00	\$3,999.00	add to cart add to list
ZLR2008	SL3 Autoclavable Handpiece	1	\$99.00	\$99.00	add to cart add to list
ZLR2004	SL3 Foot Pedal	1	\$350.00	\$350.00	add to cart add to list

Top Selling Products

New LumiSmile White
High quality whitening at an affordable price.

NEW SplashMax
Instant hydrophilicity, maximum tear strength...
Max. thickness 0.01mm.

Flashlite Magna™ 4.0
Advanced beam collimation.
Maximize curing power.

v Cart

SL3 Desktop Laser	\$3,999.00	Qty 1
SL3 Foot Pedal	\$350.00	Qty 1
SL3 Autoclavable Handpiece	\$99.00	Qty 3
Bite Registration Vanilla Bite (0-55)	\$35.75	Qty 1
Sub-Total	\$4,681.75	
Total	\$4,681.75	

checkout

v Shopping Lists

Dentist Office1
Edit add to cart

List A
Edit add to cart

LSDPfree Stuff
Edit add to cart

Your Last Complete Order
Edit add to cart

Your Last LSPD Free Order
Edit add to cart

How to Sign Up

Go to www.denmat.com. New users must select "New Customer" from the navigation bar in the upper right. To use DenMat's on-line system we must first verify your office and dental information. Please complete the form and press "Sign Me Up!"

Once the form is submitted, we will verify your information and send your login credentials to the e-mail address you provide on the form or another already on file. During working hours (5 a.m. to 5 p.m. PST) this process typically takes less than 30 minutes. On weekends and holidays, your log-in credentials will arrive early the next business day.

Please Note: Group practice or international customers, please follow your current ordering procedures, as e-commerce is not able to process these orders at this time.

Company *

License # *

First Name *

Last Name *

Address *

City *


Zip/Postal Code *

Email *

Phone *

How to Log In

Users log in by submitting a username and password in the "Log-in" field provided in the top right area of the navigation menu and pressing the "Arrow" button.

 Product Search... **Log-In**

[Forgot Password?](#) | [New Customer?](#)
1-800-4DENMAT (1-800-433-6628) | 1-805-346-3700
[Downloadable Forms](#)


[Products](#) [DenMat Lab](#) [Marketing](#) [DenMat Education](#) [Portal](#) [Fast Order](#)

How to Place an Order






DenMat recommends logging in as this allows users to see pertinent data such as "My Price, Quantity and LSDP Price," and to be able to "add to cart," "add to list," or "checkout" from the product tab or the "Fast Order" page.

There are three ways to place an order on the DenMat website.

- 1 Select a product from the top navigation bar and select "Order Now" from the product tabs and follow the checkout workflow.
- 2 Select the "Fast Order" button from the top navigation bar, choose your product tab and select "Buy Now" and follow the checkout workflow.
- 3 Use the "Product Search" tool from the top navigation bar.

 **Geristore®**

[Overview](#) [Order Now](#) [Instructions](#) [Techniques](#) [MSDS Sheets](#)
[FAQs](#) [Research](#) [Kit Options](#) [Testimonials](#)

Part#	Product Description	Qty	My Price	LSDP Price	
 030450911	Geristore A Slow Set, 10g	<input type="text" value="1"/> <input type="button" value="up"/> <input type="button" value="down"/>	\$89.00	\$66.75	<input type="button" value="add to cart"/> <input type="button" value="add to list"/>
 031457500	Geristore Syringeable Kit, Shade A2	<input type="text" value="1"/> <input type="button" value="up"/> <input type="button" value="down"/>	\$159.00	\$119.25	<input type="button" value="add to cart"/> <input type="button" value="add to list"/>
 031457520	Geristore Syringe Value Kit, Shade A3	<input type="text" value="1"/> <input type="button" value="up"/> <input type="button" value="down"/>	\$436.84	\$327.63	<input type="button" value="add to cart"/> <input type="button" value="add to list"/>
 031458550	Geristore Syringe Value Kit, Shade A2	<input type="text" value="1"/> <input type="button" value="up"/> <input type="button" value="down"/>	\$436.84	\$327.63	<input type="button" value="add to cart"/> <input type="button" value="add to list"/>
 031476110	Geristore Intro Kit	<input type="text" value="1"/> <input type="button" value="up"/> <input type="button" value="down"/>	\$77.58	\$58.19	<input type="button" value="add to cart"/> <input type="button" value="add to list"/>

Choose Shipping Method
Choose Method * 3 - 7 Days Standard \$12.95
2 - 3 Days Expedited \$17.95

Payment Options
Enter credit card information
Billing Name: ROBERT KEN DDC
Card Type: Visa
Card Number: 4111111111111111
Exp Date: 05 - 2014
CCV: 111


-- Or --
Pay On Account ☐
PO Number:

Product Search Feature

The "Product Search" feature located at the top of all pages allows users to quickly search for products by "Product Brand, Part #, Part Name or Description." Remember "logging in" offers the complete product information view including price, quantity and LSDP price.

Watch an E-Commerce Demo!

http://www.denmat.com/How_to_use_Ecommerce

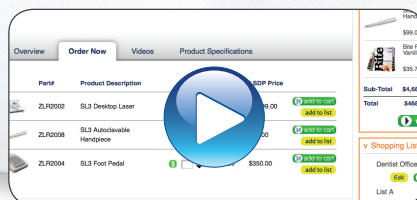

Terms & Conditions

DenMat Lab Terms & Conditions

Please see our complete terms and conditions available at www.denmat.com

Ordering

For the best possible custom-made product, please make sure you are ordering the correct size, quantity, shade, and material application. If we make an error or if your order was inaccurate, we will make every attempt to remedy the error as quickly as possible. All orders may be subject to shipping and handling charges. "Rush" orders or orders requiring special handling may be subject to additional charges.

Product Inspection

We encourage you to inspect the product prior to acceptance. You accept the product if you use it in any way, including placing the product in a patient's mouth or making modifications to it. If you do not notify the DenMat Lab of your rejection and do not return the product to the Lab within 15 business days after your receipt of the product, you will be deemed to have accepted it.

Limited Warranty

While we believe you and your patient will be delighted with our products, we understand that 100% satisfaction is not always possible. In such event, we are happy to repair or replace, at our discretion, any product that is defective or does not meet your specifications. This is your exclusive remedy. If we are at fault, there will be no charge for the repair or replacement (and no charge for shipping). If we are not at fault, there will be a 50% charge for our cost to repair or replace (in addition to the original charge for the product, plus shipping). This limited warranty is in effect from the date you received the product and for a period of 90 days from that date.

For all DenMat Lab restorations (excluding Snap-On Smile), DenMat's warranty obligation, with respect to an original restoration affixed to a particular tooth, is limited to one (1) time replacement for that restoration within five (5) years.

For Snap-On Smile, the warranty obligation is limited to one (1) time replacement of the original Snap-On Smile (both full and partial arches) within one (1) year and covers any defects in materials or workmanship.

These limited warranties are in effect from the date you received the product and are between DenMat and the doctor. There will be no handling charge for warranty work for the first 90 days of the warranty period. After the 90 days, there will be a \$37 handling fee (US dollars) plus shipping (depending on current shipping rates) on all warranty claims.

This limited warranty is in lieu of all other warranties, express or implied and is void if the product is improperly stored or used. There are no implied warranties of merchantability, fitness for a particular purpose or otherwise. Before using this product, you agree that it is your responsibility to determine whether it is suitable for the intended use and will assume all risk and/or liability associated with use of or inability to use the product. This warranty gives you specific legal rights, and you may also have other rights that vary from state to state. Some states do not allow the exclusion or limitation of incidental or consequential damages so this limitation or exclusion may not apply to you.

DenMat Professional, Consumable & Equipment Terms & Conditions

Please see our complete terms and conditions available at www.denmat.com

Ordering

So that we may provide you the best service possible, please make sure you are ordering the correct size, quantity, shade, and material application. If we make an error or if your order was inaccurate, we will make every attempt to remedy the error as quickly as possible. All orders may be subject to shipping and handling. "Rush" orders or orders requiring special handling may be subject to additional charges.

Cancellations

To cancel an order that has not shipped, please call your Customer Service Representative. If your order has not shipped, we will cancel it. If your order has shipped, please obtain a Return Authorization Number from your Customer Service Representative because we cannot accept returns without prior authorization. Return Shipping Labels will be provided to you if an order was sent to you in error.

Limited Warranty

Professional Consumable Products

While we believe you and your patients will be delighted with our products, we understand that 100% satisfaction is not always possible. In such event, we are happy to replace or refund your purchase price of any DenMat professional consumable product that is defective if you notify us within 30 days of the purchase date. This is your exclusive remedy.

Equipment

We are happy to diagnose and repair DenMat equipment found defective within the equipment warranty period. Equipment outside of the warranty period will be evaluated and subject to the results of evaluation. DenMat will contact you with a price quotation for the cost of repairs or replacement, if applicable. Any repairs may be made with new or refurbished parts at the sole discretion of DenMat.

- Shipping and handling charges for equipment within warranty will be paid for by DenMat.
- Shipping and handling charges for equipment outside warranty and for professional consumable products is the responsibility of the customer.

This is a limited warranty (for both Professional consumable Products and Equipment) and is in lieu of all other warranties, expressed or implied, and is void if the product is improperly stored or used. There are no implied warranties of merchantability, fitness for a particular purpose or otherwise. Before using this product, you agree that it is your responsibility for determining whether it is suitable for the intended use and will assume risk and liability associated with use of the product. You agree that DenMat will not be liable for any economic, incidental or consequential loss or damage that arises out of the use of or the inability to use its products. This warranty gives you specific legal rights, and you may also have other rights that vary from state to state. Some states do not allow the exclusion or limitation of incidental or consequential damages so this limitation or exclusion may not apply to you.

Payment Terms

We appreciate payment on your account within 30 days of Invoice date. In the event your account balance is past due, we reserve the right to impose a service charge of 1.5% per month to your account until payment is brought current. In the event your account balance is substantially past due and we have no choice but to refer your account to another party for collection, please be aware that you will be responsible for any collection fees and expenses.

Product Return Policy

Consumables Returns (porcelain, materials, etc.):

- Product returned within 30 days of purchase will be accepted for full refund. No questions asked.
- If the product is defective and you notify us after 30 days, we would be happy to credit your account upon our confirmation that the product is defective.
- We are unable to accept returns for products beyond 60 days. Returns received 31– 60 days will be given a 50% credit return.
- Please be aware that your account may be subject to shipping charges.
- Promotional coupon credit has an expiration date of six (6) months from time of redemption.
- Concession credits have an expiration date of six (6) months from time of application.
- Incorrect product shipped to a customer: Products will be promptly replaced with correct product upon DenMat issuing RMA number. If the incorrect products have not been received by DenMat within 30 days of issuance of the RMA number, the customer will be charged for the products. All returned products must be in original, unopened condition.
- Damaged products: Products that have been damaged in-transit to the customer will be replaced, within 30 days of receipt by the customer, at no charge to the customer by DenMat.

Equipment Returns:

- Product returned within 30 days of purchase will be accepted for full refund. No questions asked.
- We are happy to repair or replace equipment within its warranty period. DenMat retains sole discretion to determine whether repair or replacement is appropriate.
- Equipment within warranty will be repaired according to the terms of the specific equipment warranty.
- Equipment outside warranty will be repaired according to the terms of the specific equipment warranty, with all charges (including parts, labor) paid by the customer.
- Please be aware that your account may be subject to shipping charges.
- Equipment returned for credit 31 – 60 days from date of purchase will be subject to a 15% restocking fee.
- Equipment returned within 30 days from date of purchase will be given full credit of the purchase price.
- Incorrect product shipped to a customer: Products will be promptly replaced with correct product upon DenMat issuing RMA number. If the incorrect products have not been received by DenMat within 30 days of issuance of the RMA number, the customer will be charged for the products. All returned products must be in original, unopened condition.
- Damaged products: Products that have been damaged in-transit to the customer will be replaced, within 30 days of receipt by the customer, at no charge to the customer by DenMat.

Indemnity

You agree to hold DenMat and its respective officers, directors, employees, agents, successors, and assigns harmless from and against any claim, loss, damage or expense, including attorney's fees and court costs, arising from or related to your use or misuse of DenMat's products.

Our Contact Information

If you have questions or comments about these Terms and Conditions, we are available Monday – Friday from 6 a.m. – 5 p.m. (Pacific Time). Please feel free to contact our Customer Service team in any of three ways:

1. Send us an e-mail at customerservice@denmat.com
2. Call us at 1-800-433-6628
3. Write us at:
Den-Mat Holdings, LLC
Attn: Customer Service
P.O. Box 1729
Lompoc, CA 93438-1729

Contact Us

Phone

1-800-4DENMAT (1-800-433-6628)
1-805-346-3700

Street

Den-Mat Holdings, LLC
1017 W. Central Ave.
Lompoc, CA 93436-2701

Mail

Den-Mat Holdings, LLC
P.O. Box 1729
Lompoc, CA 93438-1729

E-mail

Please e-mail your comments,
suggestions or questions to
info@denmat.com

Web

Visit us on-line at www.denmat.com



©2016 Den-Mat Holdings, LLC. All rights reserved. ATRIDOX is a registered trademark of Tolmar, Inc. Adidas is a registered trademark of Adidas AG. Hogies is a registered trademark of Martin Hogan Pty Ltd. VELscope VX is a registered trademark of LED Medical Diagnostics, Inc. IPS e.Max and Zirpress are registered trademarks of Ivoclar Vivadent, Inc. Captek is a registered trademark of Precious Chemicals Ltd. Aquasil is a trademark of Dentsply International, Inc. Atlantis is a registered trademark of Dentsply IH Inc. Imprint is a registered trademark of 3M Company. Bellatek and Encode are registered trademarks of Biomet 3i, LLC. ARESTIN is a registered trademark of Orapharm, Inc. Orascope is a registered trademark of Kerr Corporation. Surgitel is a registered trademark of General Scientific Corporation. Zeiss is a registered trademark of Carl Zeiss Ah. Rembrandt is a registered trademark of Johnson and Johnson Corporation.
820048630 07/16SN



1017 W. Central Avenue, Lompoc, CA 93436
1-800-4DENMAT (1-800-433-6628)

