





Snap-On Smile redefines provisionals. It's non-invasive, it's easy and your patients will love it. No more messy temporaries or complicated procedures requiring too much of your chair time. Snap-On Smile's range of **versatile clinical applications** includes **Opening & Restoring Vertical Dimension, Removable Partial Denture, Implant Provisional and Attachment, and Esthetic/Cosmetic Smile Enhancement** (including alignment, whitening and severe crowding). Snap-On Smile is available as both a full arch (6 units or more) and as a partial arch (5 units or less).



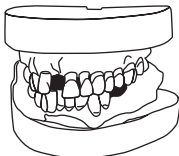
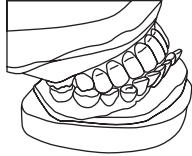

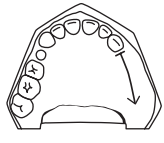

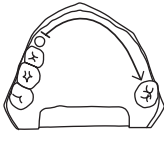
SNAP-ON SMILE INDICATIONS & CONTRAINDICATIONS

INDICATIONS

BEFORE & AFTER IMAGE

Opening & Restoring Vertical Dimension	 <p>Patient: Jenny G.</p>
Removable Partial Denture (Missing Teeth)	 <p>Patient: Yolanda S.</p>
Implant Provisional and Attachment (Embedding)	 <p>Patient: Grizelda G.</p>
Esthetic/Cosmetic Smile Enhancement (Alignment, Whitening, Severe Crowding)	 <p>Patient: Sharon B.</p>



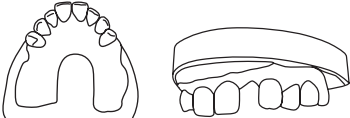

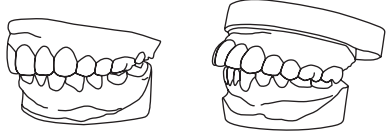
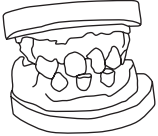

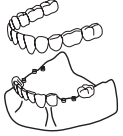

CONTRAINDICATIONS

Severe Periodontal Cases (Generalized Mobile Teeth)		Severe Class III Bites	
Edentulous		Cantilever Distal Extensions > 22mm	
Snap-On Smile Over an Existing Denture		Edentulous Spans > 40mm	

SNAP-ON SMILE STANDARD VS. ADVANCED CASES – PRODUCT PARAMETERS

Standard Snap-On Smile cases include those noted below. For more difficult, complex cases, **Snap-On Smile Advanced** is the optimal alternative, for which **Digital Preview Technology** and **Premium Consult** services allow you to ensure that the design of your more severe cases are made to your precise criterion. Snap-On Smile Advanced takes the guesswork out of complex case design and offers **predictability** every time.

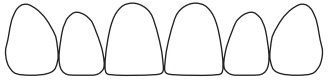





<ul style="list-style-type: none"> Edentulous spans up to 22mm (Full Arch) & up to 14mm (Partial Arch) 	 <p>Full Arch Partial Arch</p>	<ul style="list-style-type: none"> Edentulous spans with 2 abutments > 22mm & < 40mm (Full Arch) 	
<ul style="list-style-type: none"> Up to 22mm cantilever extensions (Full Arch) & up to 9mm cantilever extensions (Partial Arch) 	 <p>Full Arch Partial Arch</p>	<ul style="list-style-type: none"> Correcting severe buccal inclinations 	
<ul style="list-style-type: none"> Appropriate for Class I (Ideal Bite) or Class II (Overbite) 	 <p>Class I (Ideal Bite) Class II (Overbite)</p>	<ul style="list-style-type: none"> Severe malocclusion (Cross bite) 	
<ul style="list-style-type: none"> Appropriate for certain Class III bites (Underbites) 		<ul style="list-style-type: none"> Embedded Implants (Up to 3) 	
<ul style="list-style-type: none"> Cases with limited retention 			

SNAP-ON SMILE SMILE STYLES

Snap-On Smile cases can be designed in 5 customized Smile Styles based on your preference for the shape of the Centrals, Laterals and Cuspids.

SNAP-ON SMILE SMILE STYLES

<ul style="list-style-type: none"> Hollywood* (<u>Centrals</u>: Square Round, <u>Laterals</u>: Square Round, <u>Cuspids</u>: Flat (with option to make slightly shorter)) 	
<ul style="list-style-type: none"> Softened* (<u>Centrals</u>: Square Round, <u>Laterals</u>: Square Round, <u>Cuspids</u>: Round) 	
<ul style="list-style-type: none"> Natural* (<u>Centrals</u>: Square Round, <u>Laterals</u>: Short Square Round, <u>Cuspids</u>: Pointed) 	
<ul style="list-style-type: none"> Functional* (<u>Centrals</u>: Square Round, <u>Laterals</u>: Square Round, <u>Cuspids</u>: Pointed) 	
<ul style="list-style-type: none"> Follow Existing Dentition (Default) 	